

Catecholamines in 100% Aqueous by Gemini 5um C18

Page 1 of 1

Column: Gemini 5u C18 110A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic


Eluent A: 20mM Phosphate buffer pH2.5

Pct. A: 100


Flow Rate: 1.5 mL/min

Col. Temp.: 30 °C

Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

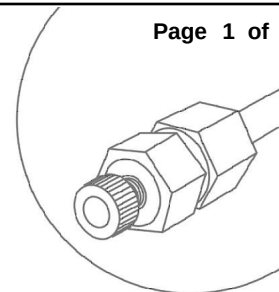

Products used in this Application:

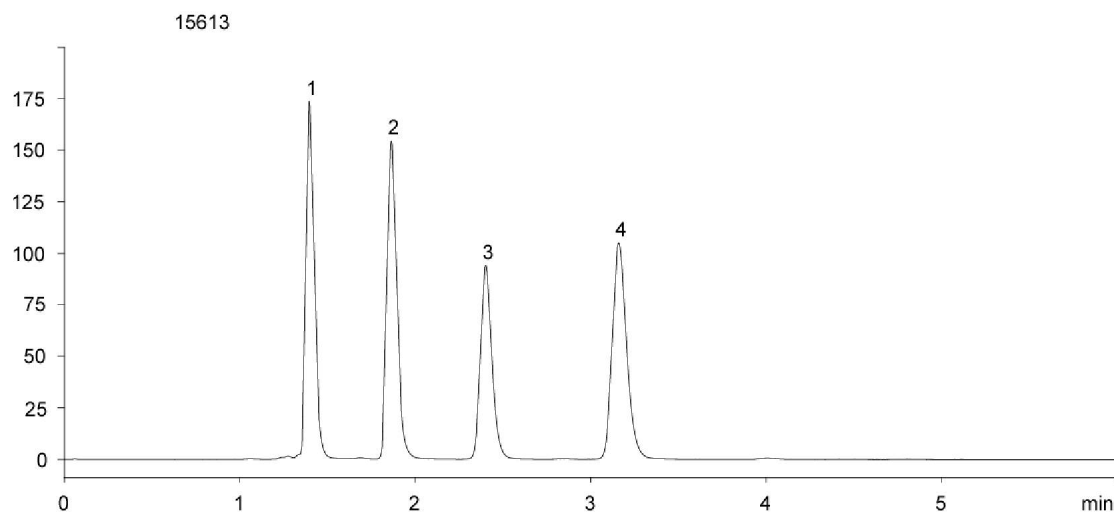

ANALYTES:

- 1 Norepinephrine
- 2 Epinephrine
- 3 Normetanephine
- 4 Dopamine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.