

Food Coloring

Column: Luna 5u C18(2) 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Isocratic

Eluent A: 5 mM TBA Phosphate buffer pH 7.0

Pct. A: 55

Eluent B: Methanol


Pct. B: 45

Flow Rate: 0.5 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 630/610/502/428/484 nm (ambient)

Page 1 of 1

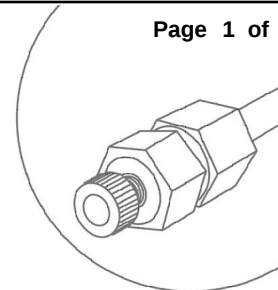

Products used in this Application:

LUNA™

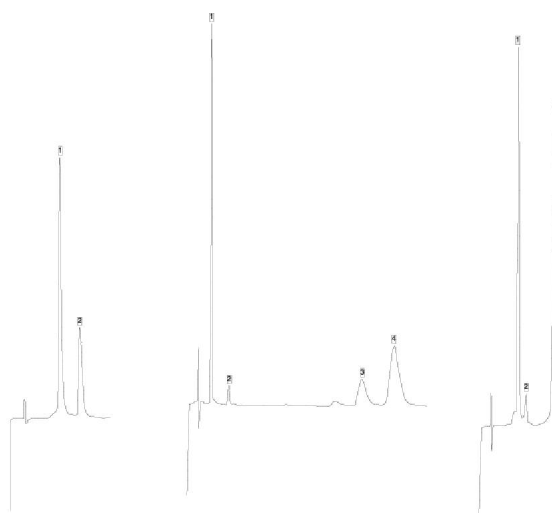

ANALYTES:

- 1 Yellow 5, FDC
- 2 Yellow 6, FDC
- 3 Blue 2, FDC
- 4 Impurity
- 5 Impurity
- 6 Blue 1, FDC
- 7 Red 2, FDC
- 8 Red 4, FDC
- 9 Red 40, FDC



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.