


Triazine Herbicides: EPA Method 619

Column: Zebtron ZB-5ms
Dimensions: 30meter x 0.25mm x 0.25µm
Order No.: 7HG-G010-11
Injection: Split 20:1 @ 250°C, 1µL
Carrier Gas: Constant Flow Helium 1.2 mL/min
Oven Program: 140C to 220°C at 4°C/min
Detector: MSD; 30-350amu
Sample:

Gold Seal: Phenomenex single groove (AGO-7518)
Liner: Single Taper Gooseneck (AGO-4680)
Septum: Phenomenex PhenoRed-400 (AGO-4690)

1. Atraton
 2. Simazine
 3. Prometon
 4. Atrazine
 5. Propazine
 6. Tertbutylazine
 7. Secbumeton
 8. Simetryn
 9. Ametryn
 10. Prometryn
 11. Terbutryn

Notes: Analytes are 100ppm.

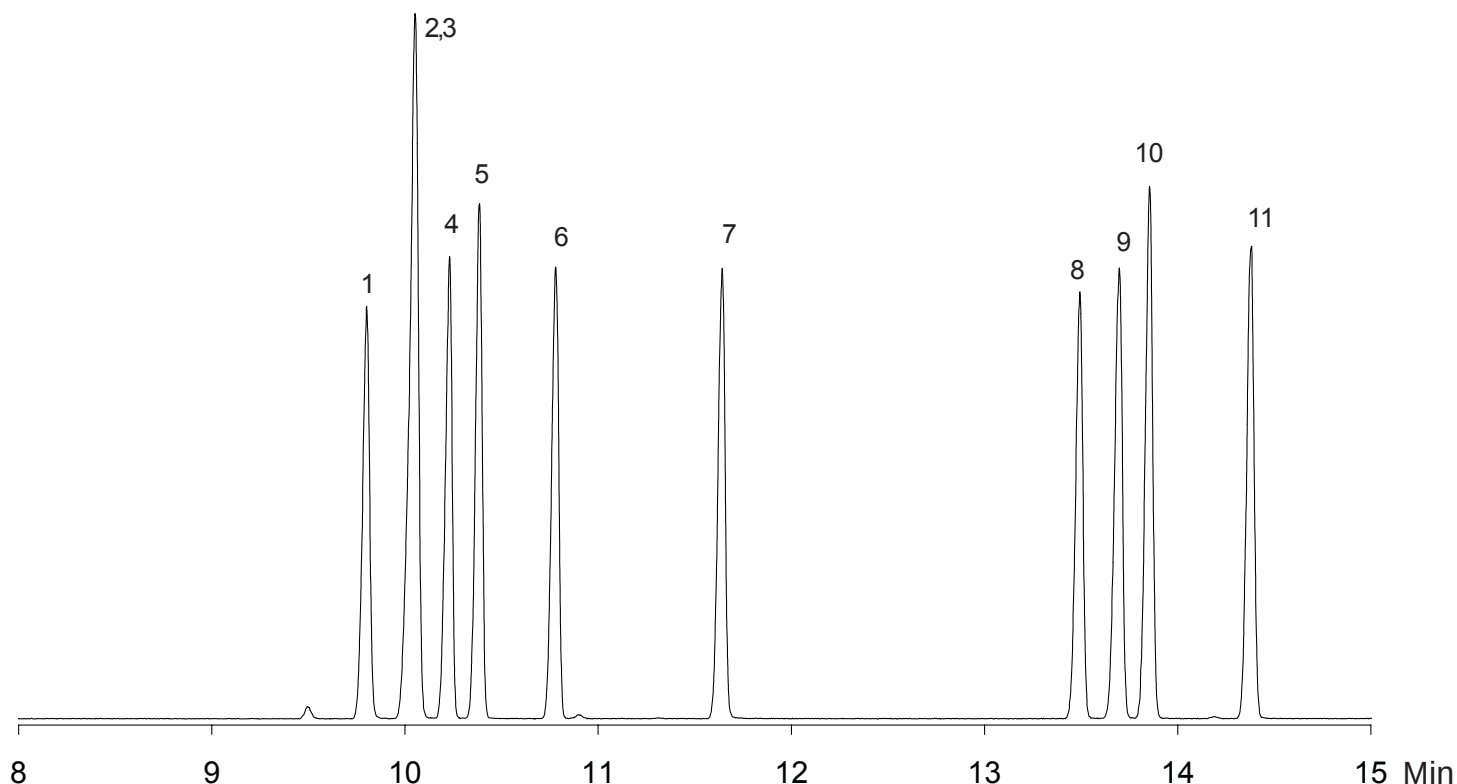

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.



 ...breaking with traditionSM

USA	Puerto Rico	Canada	United Kingdom	Germany	New Zealand	Australia
mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA	271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926	411 Madrid Ave. Torrance, CA 90501-1430 USA	Queens Avenue, Hurdsfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England	Zeppelinstr. 5 63741 Aschaffenburg Deutschland	P.O. Box 31-601 Milford Auckland New Zealand	P.O. Box 264 Pennant Hills NSW 1715 Australia
tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	(800) 541-HPLC (310) 328-7768 info@phenomenex.com	(800) 543-3681 (310) 328-7768 info@phenomenex.com	01625-501367 01625-501796 ukinfo@phenomenex.com	06021-58830-0 06021-58830-11 anfrage@phenomenex.com	09-4780951 09-4780952 info@phenomenex.co.nz	1800-553-929 1800-553-923 info@phenomenex.com.au