


Onyx RP (C18) Test Chromatogram (AL0-7836)

Column: Onyx Monolithic C18**Dimensions:** 50 x 4.6 mm ID**Order No:** CH0-7644**Elution Type:** Isocratic**Eluent A:** Acetonitrile**Pct. A:** 60**Eluent B:** Water**Pct. B:** 40**Flow Rate:** 2 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)


Page 1 of 1

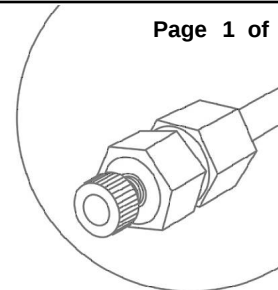

Products used in this Application:

ONYX

Analyst Note: This standard is used to test the following Phenomenex Onyx™ Monolithic Reversed Phase Columns: Onyx C8 and C18.
NOTE: Separation conditions and results may vary for your column. Contact your Technical Representative.
Efficiency (based on Anthracene) = 92136 p/m (Test Mix Order No.: 7836)

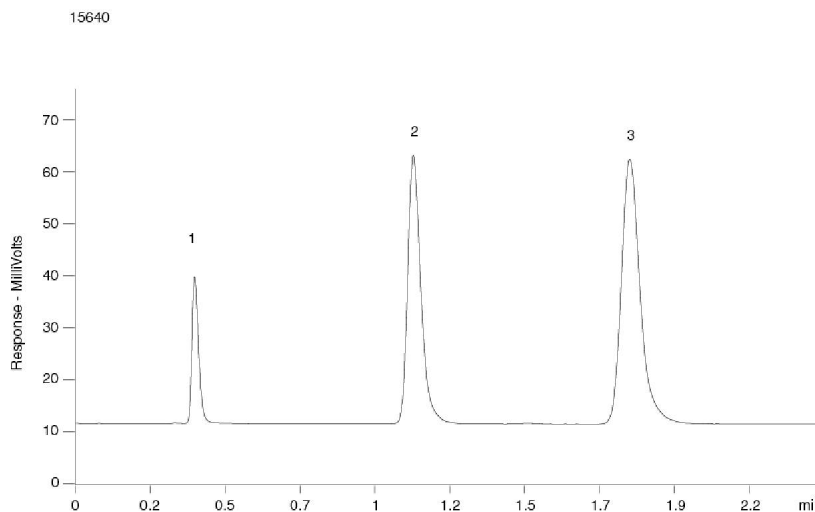**ANALYTES:**

- 1 Thiourea
- 2 Progesterone
- 3 Anthracene

AUTHORIZED DISTRIBUTOR

SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr