


Organochlorine Pesticides: EPA Method 8081

Column:	Zebron ZB-5ms	Gold Seal:	Phenomenex single groove	(AG0-7518)
Dimensions:	30meter x 0.25mm x 0.25µm	Liner:	Split/Splitless Single Taper	(AG0-7516)
Order No.:	7HG-G010-11	Septum:	Phenomenex PhenoRed-400	(AG0-4690)
Injection:	Split 20:1 @ 250°C, 1µL			
Carrier Gas:	Constant Flow Helium 1.0 mL/min			
Oven Program:	100°C to 200°C at 20°C/min to 220°C at 5°C/min to 300°C at 15°C/min for 2 min			
Detector:	MSD; 30-400 amu			
Sample:	1. Etridiazole 2. Chloroneb 3. Propachlor 4. Trifluralin 5. α-BHC (Lindane) 6. Hexachlorobenzene 7. β-BHC 8. γ-BHC 9. Chlorothalonil 10. δ-BHC 11. Heptachlor 12. Chlorpyrifos 13. DCPA (Dacthal) 14. Aldrin 15. Heptachlor Epoxide 16. γ-Chlordane 17. α-Chlordane 18. Endosulfan I 19. 4,4'-DDE 20. Dieldrin 21. Endrin 22. Chlorobenzilate 23. Endosulfan II 24. 4,4'-DDD 25. Endrin Aldehyde 26. Endosulfan Sulfate 27. 4,4'-DDT 28. Endrin Ketone 29. Methoxychlor 30. cis-Permethrin 31. trans-Permethrin			

Notes: All analytes are 100 ppm with the exception of 30 and 31 at a combined concentration of 200 ppm and 16 and 17 at 200 ppm each.

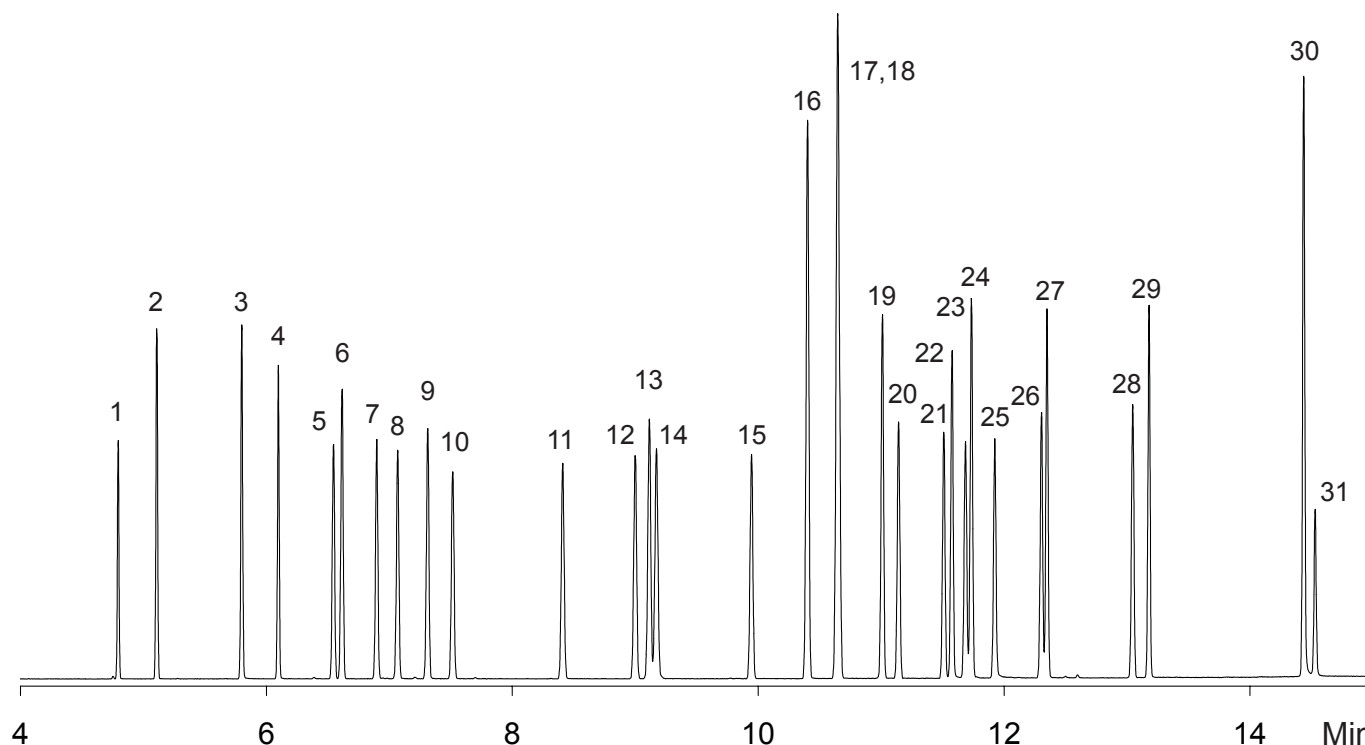

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex
...breaking with tradition™

USA
mail: 411 Madrid Avenue
Torrance, CA
90501-1430
USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
273 Sierra Morena,
Suite #104
San Juan,
Puerto Rico 00926
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
411 Madrid Avenue
Torrance, CA
90501-1430
USA
(800) 543-3681
(310) 328-7768
info@phenomenex.com

United Kingdom
Queens Avenue,
Hurdfield Ind. Est.,
Macclesfield, Cheshire
SK10 2BN, England
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
Zeppelinstr. 5
63741 Aschaffenburg
Deutschland
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
P.O. Box 31-601
Milford
Auckland
New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
P.O. Box 4084
Lane Cove, NSW 1595
Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au