


APPLICATION NOTE - HPLC
Amoxicillin - Alternative Method

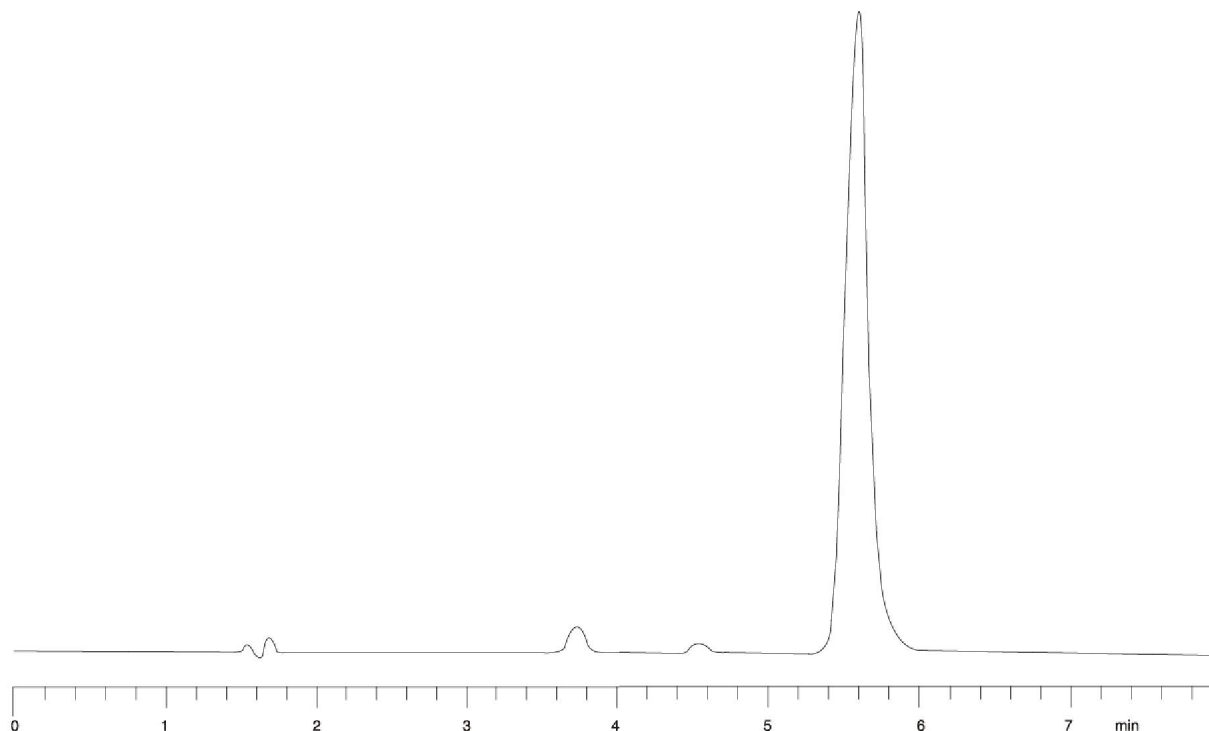

Column: Luna 5 μ m C8 (2)

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4249-E0

Elution Type: Isocratic

Eluent A: 20mM Ammonium acetate, pH 4.9

Pct. A: 95

Eluent B: Methanol

Pct. B: 5

Flow Rate: 1 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 230 nm

SAMPLE ANALYTES:

1 Amoxicillin

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 3294/USP/72

Page 1 of 1