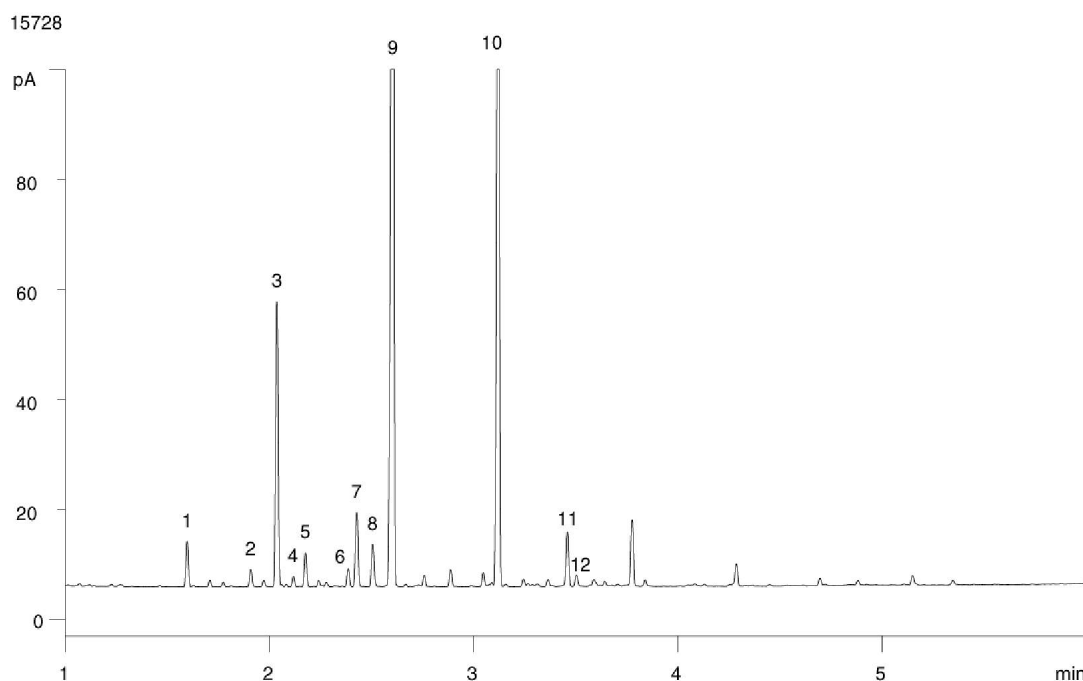


Amino Acid Analysis of Apple Juice by EZ:faast

Page 1 of 1

Column: EZ:faast GC
Phase: Zebron ZB-AAA (proprietary)
Dimensions: 10 meters x 0.25 mm x 0 µm
Order No: KG0-7165
Oven Profile: 110 to 320°C at 32°C/min
Carrier Gas: Constant Flow Helium, 1.5 mL/min
Injection: Split 15:1 2 µL @ 250 °C
Detection: Flame Ionization (FID) (320 °C)

Products used in this Application:

EZ:faast™

ANALYTES:

- | | |
|-------------|--------|
| 1 ALA | 10 ASP |
| 2 VAL | 11 GLU |
| 3 NORVALINE | 12 PHE |
| 4 LEU | |
| 5 ILE | |
| 6 THR | |
| 7 SER | |
| 8 PRO | |
| 9 ASN | |

©Copyright 2000-2005 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with tradition™

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au