

Deoxynucleosides using Synergi Fusion-RP

Page 1 of 1

Column: Synergi 4u Fusion-RP 80A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4424-E0

Elution Type: Gradient

Eluent A: 20mM Ammonium acetate, pH 4.5

Eluent B: Acetonitrile


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	10	60	40

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 260 nm (ambient)

Products used in this Application:

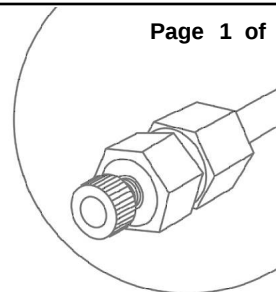synergi[™]

15766

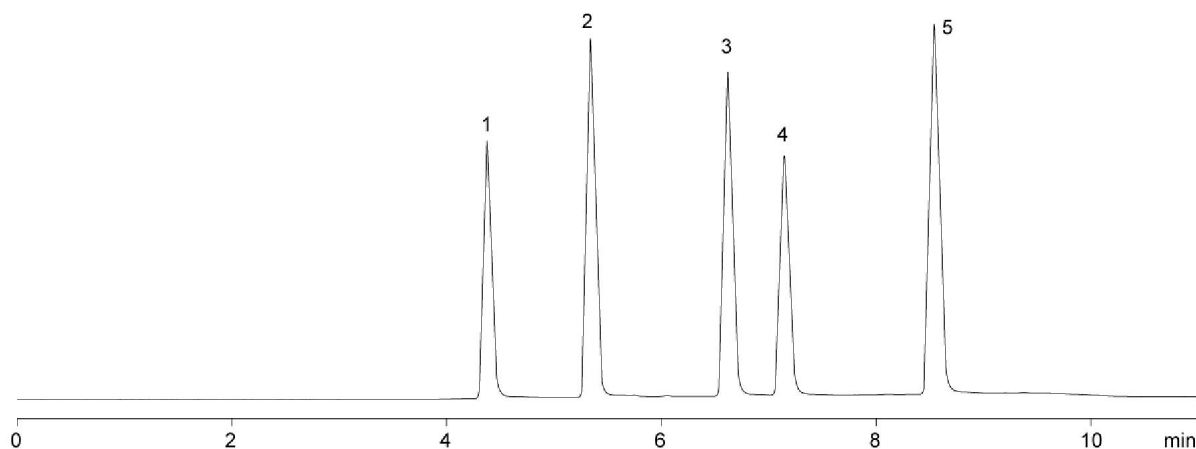

ANALYTES:

- 1 2'-Deoxycytidine
- 2 2'-Deoxyuridine
- 3 2'-Deoxyguanosine
- 4 2'-Deoxythymine
- 5 2'-Deoxyadenine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.