


Kerosene

Column: Zebtron ZB-5
Dimensions: 30m x 0.53mm x 1.5µm
Order No.: 7HK-G002-28
Injection: Split 30:1 @ 275°C, 1µL
Carrier Gas: Constant Flow Helium 4.1mL/min
Oven Program: 40°C for 5 min to 300°C at 4°C/min for 5 min.
Detector: FID @ 300°C
Sample: 1. Kerosene

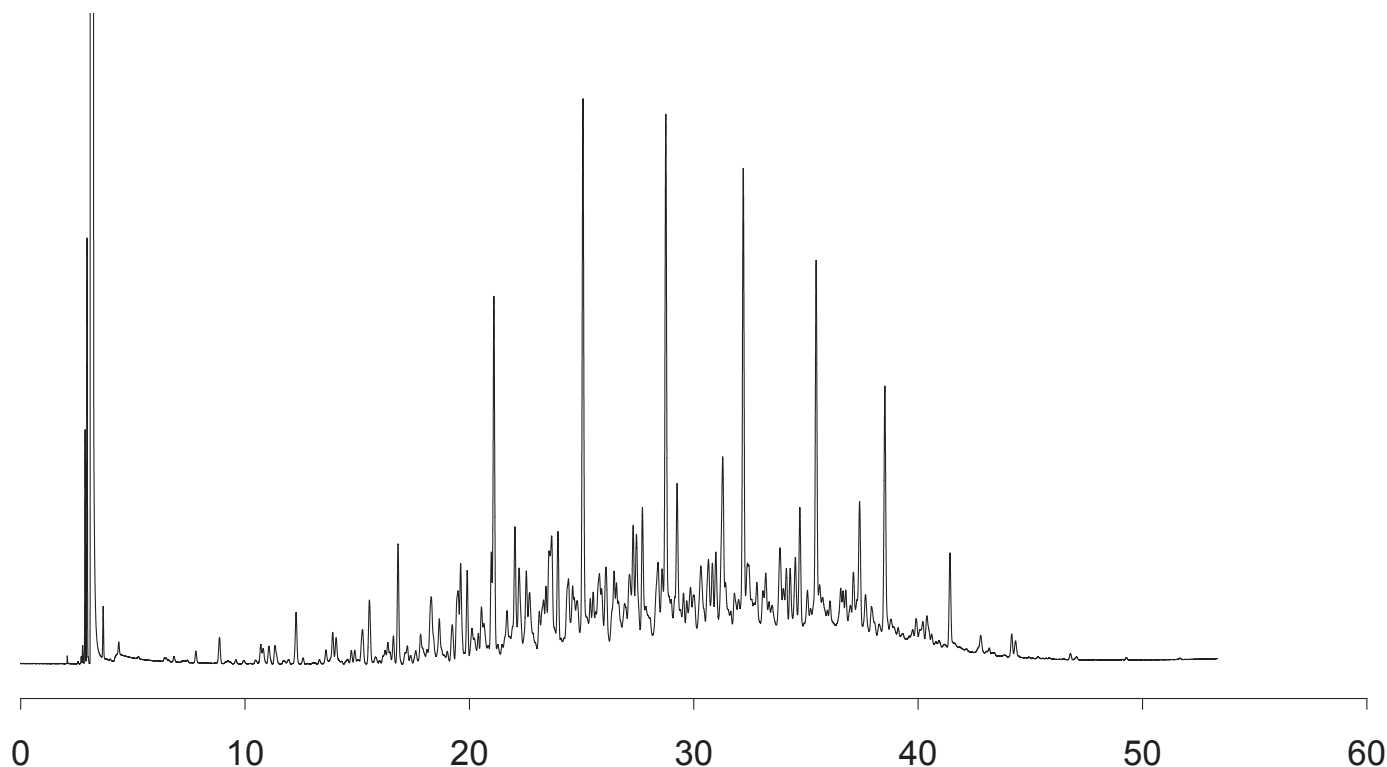

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

 **phenomenex**[®]
*...breaking with tradition*SM

USA	Puerto Rico	Canada	United Kingdom	Germany	New Zealand	Australia
mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA	271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926	411 Madrid Ave. Torrance, CA 90501-1430 USA	Queens Avenue, Hurdsfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England	Zeppelinstr. 5 63741 Aschaffenburg Deutschland	P.O. Box 31-601 Milford Auckland New Zealand	P.O. Box 264 Pennant Hills NSW 1715 Australia
tel.: (310) 212-0555	(800) 541-HPLC	(800) 543-3681	01625-501367	06021-58830-0	09-4780951	1800-553-929
fax: (310) 328-7768	(310) 328-7768	(310) 328-7768	01625-501796	06021-58830-11	09-4780952	1800-553-923
email: info@phenomenex.com	info@phenomenex.com	info@phenomenex.com	ukinfo@phenomenex.com	anfrage@phenomenex.com	info@phenomenex.co.nz	info@phenomenex.com.au