


Phthalate Esters - EPA Method 606

Column: Zebron ZB-5ms
Dimensions: 30meter x 0.32mm x 0.25µm
Order No.: 7HM-G010-11
Injection: Split 20:1 @ 300°C, 1µL
Carrier Gas: Constant Flow Helium 2.6mL/min
Oven Program: 75°C to 300 at 15°C/min for 5 min
Detector: FID @ 300°C
Sample:

1. Dimethyl Phthalate
2. Diethyl Phthalate
3. Di-n-butyl Phthalate
4. Butyl Benzyl Phthalate
5. Bis(2-ethylhexyl)phthalate
6. Di-n-octyl Phthalate

Gold Seal: Phenomenex single groove (AGO-7518)
Liner: Single Taper Gooseneck (AGO-4680)
Septum: Phenomenex PhenoRed-400 (AGO-4690)

Notes: Analytes are 100ppm in methanol.

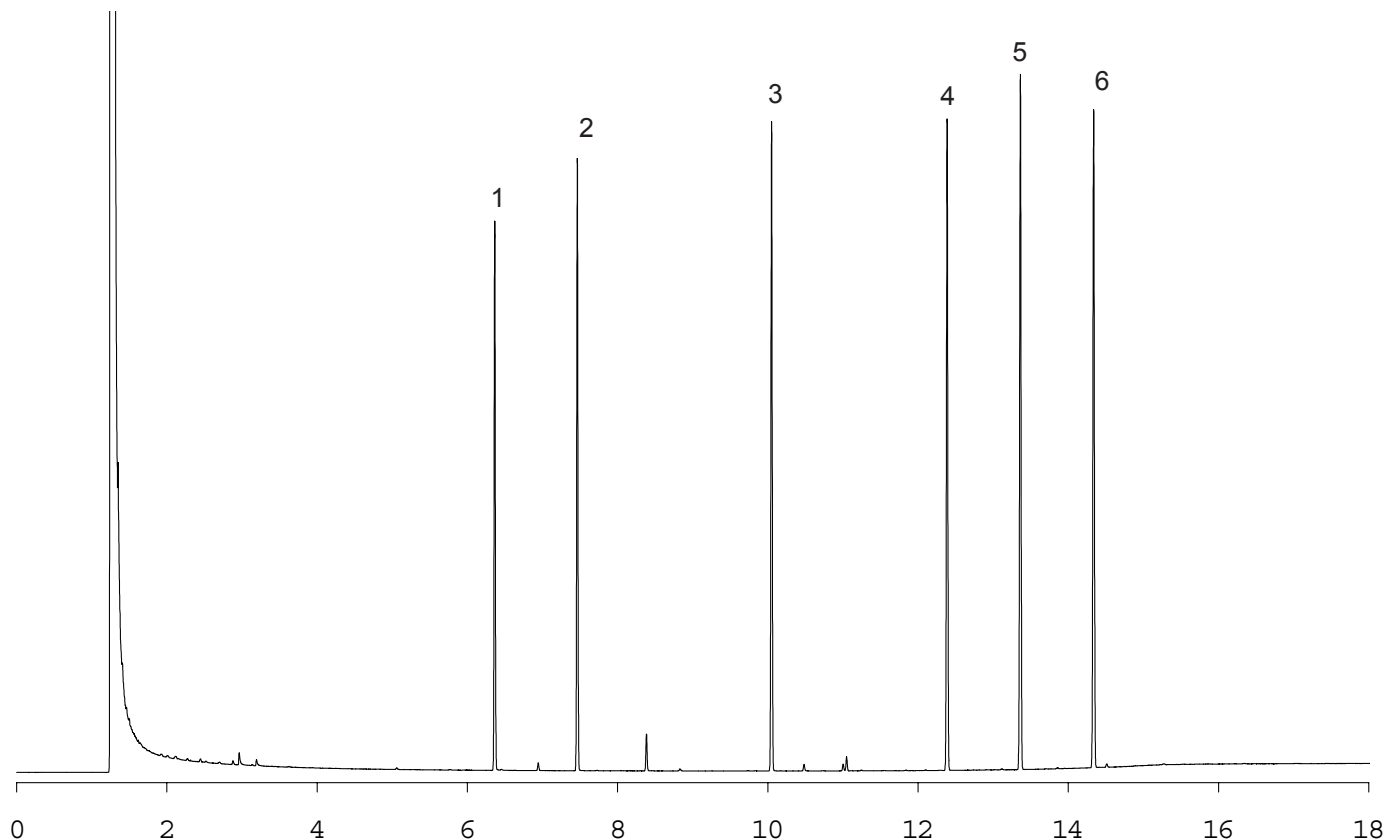

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.



...breaking with traditionSM

USA mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	Puerto Rico 271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926 (800) 541-HPLC (310) 328-7768 info@phenomenex.com	Canada 411 Madrid Ave. Torrance, CA 90501-1430 USA (800) 543-3681 (310) 328-7768 info@phenomenex.com	United Kingdom Queens Avenue, Hurdsfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England 01625-501367 01625-501796 ukinfo@phenomenex.com	Germany Zeppelinstr. 5 63741 Aschaffenburg Deutschland 06021-58830-0 06021-58830-11 anfrage@phenomenex.com	New Zealand P.O. Box 31-601 Milford Auckland New Zealand 09-4780951 09-4780952 info@phenomenex.co.nz	Australia P.O. Box 264 Pennant Hills NSW 1715 Australia 1800-553-929 1800-553-923 info@phenomenex.com.au
--	--	---	---	---	---	--