

Immunoglobulins on BioSep-SEC-S3000

Column: BioSep-SEC-S 3000

Dimensions: 300 x 7.8 mm ID

Order No: 00H-2146-K0

Elution Type: Isocratic


Eluent A: 50 mM NaH₂PO₄, pH 6.8

Pct. A: 100

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 280 nm (ambient)

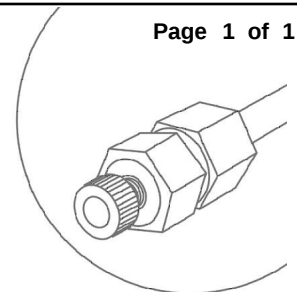

Products used in this Application:

BioSep™

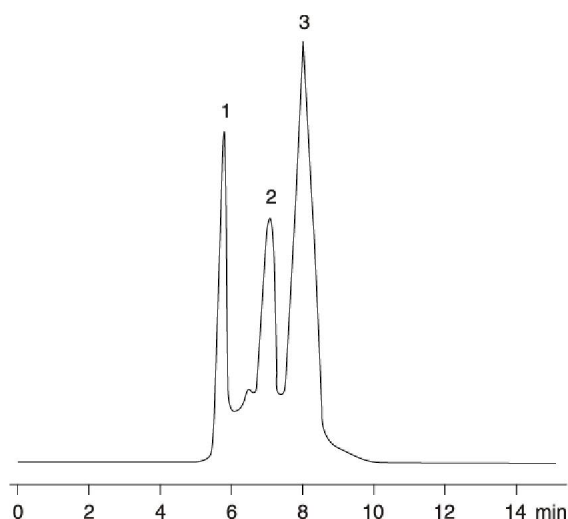

ANALYTES:

- 1 IgM
- 2 IgA
- 3 IgG



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.