


USP 467 Residual Solvents

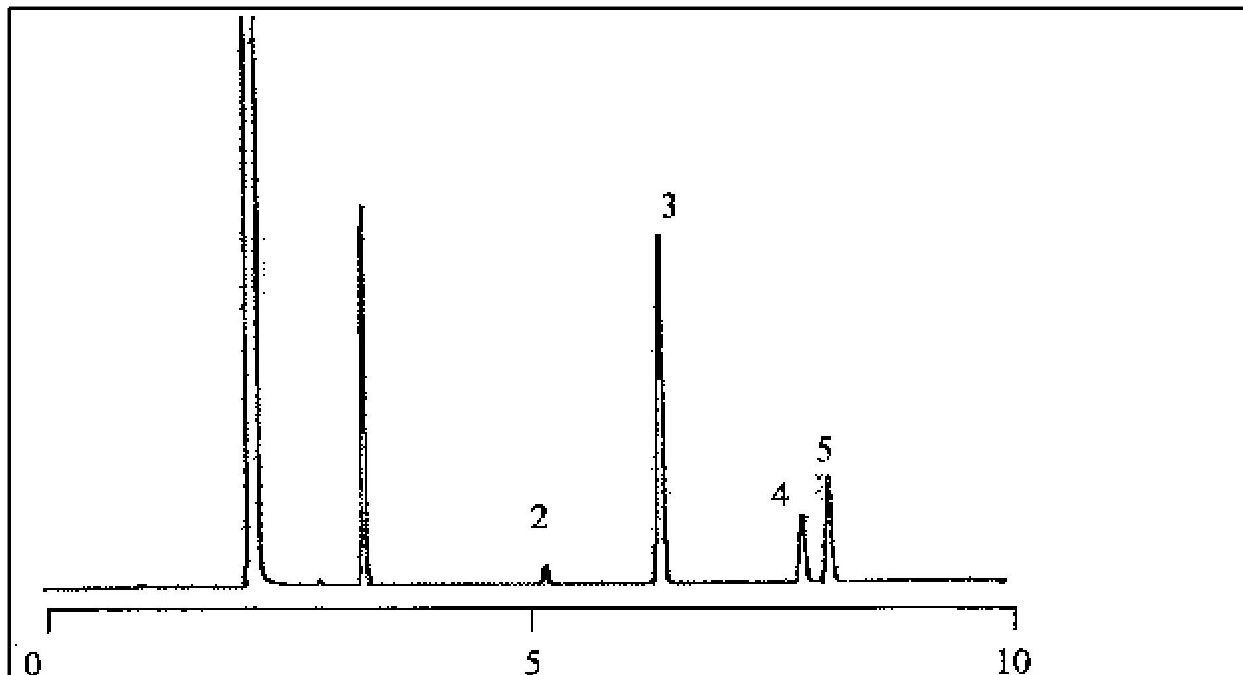

Column: Zebron ZB-5
Phase: 5%-diphenyl-95%-dimethyl polysiloxane copolymer
Dimensions: 30 meters x 0.53 mm x 1.5 μ m
Order No:

Gradient Profile:

Step No.	Step Time	Flow Rate	Temperature
1	0	0	40
2	4	0	40
3	14	0	80

Carrier Gas: Helium 25 cm/sec @ 40 °C
Flow Rate: 3.31 mL/min
Injection: Split 100:1 1 μ L @ 250 °C
Detection: Flame Ionization (FID) 300 °C

SAMPLE ANALYTE

- 1 Methylene chloride
- 2 Chloroform
- 3 Benzene
- 4 Trichloroethylene
- 5 1,4-Dioxane

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. 1267/ZEB/007
No.