


Nitrogen Pesticides - EPA Method 633

Column:	Zebtron ZB-5ms	Gold Seal:	Phenomenex single groove	(AG0-7518)
Dimensions:	30meter x 0.25mm x 0.25µm	Liner:	Split/Splitless Single Taper	(AG0-7516)
Order No.:	7HG-G010-11	Septum:	Phenomenex PhenoRed-400	(AG0-4690)
Injection:	Split 20:1 @ 250°C, 1µL			
Carrier Gas:	Constant Flow Helium 1.0mL/min			
Oven Program:	100°C to 200°C at 20°C/min to 220°C at 5°C/min to 300°C at 15°C/min for 2 min			
Detector:	MSD; 30-400amu			
Sample:	1. Deet 2. Terbacil 3. Metribuzin 4. Bromacil 5. Triadimefon (amiral) 6. Tricyclazole 7. Hexazinone			

Notes: Analytes are 100ppm.

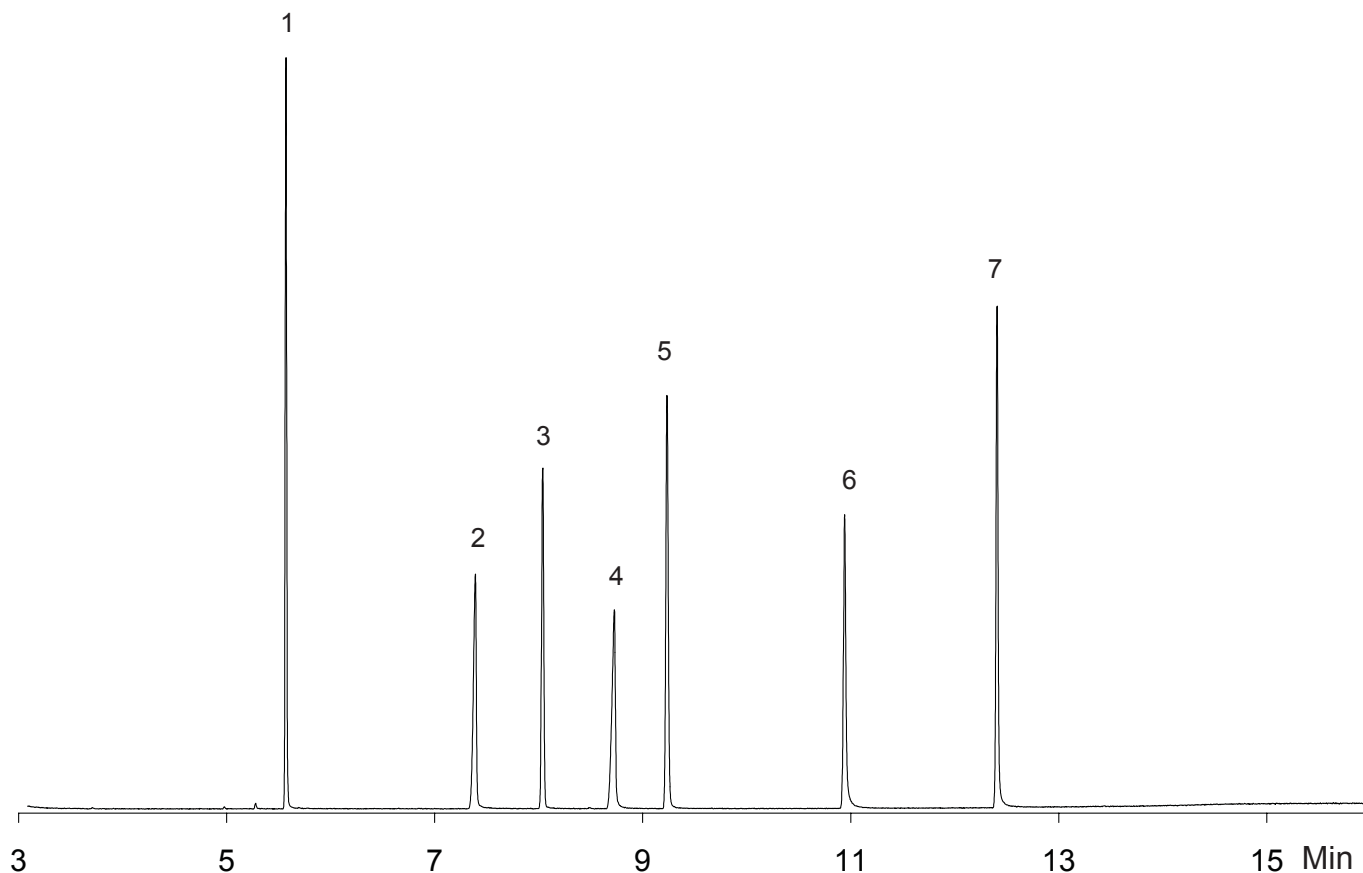

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex
...breaking with tradition™

USA	Puerto Rico	Canada	United Kingdom	Germany	New Zealand	Australia
mail: 411 Madrid Avenue Torrance, CA 90501-1430 USA	273 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926	411 Madrid Avenue Torrance, CA 90501-1430 USA	Queens Avenue, Hurdfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England	Zeppelinstr. 5 63741 Aschaffenburg Deutschland	P.O. Box 31-601 Milford Auckland New Zealand	P.O. Box 4084 Lane Cove, NSW 1595 Australia
tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	(800) 541-HPLC (310) 328-7768 info@phenomenex.com	(800) 543-3681 (310) 328-7768 info@phenomenex.com	01625-501367 01625-501796 ukinfo@phenomenex.com	06021-58830-0 06021-58830-11 anfrage@phenomenex.com	09-4780951 09-4780952 info@phenomenex.co.nz	02-9428-6444 02-9428-6445 info@phenomenex.com.au