

Proteins

Column: Jupiter 5u C18 300A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4053-E0

Elution Type: Gradient

Eluent A: 0.1% TFA in water


Eluent B: 0.1% TFA in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25
	2	15	45	55


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 220 nm (ambient)

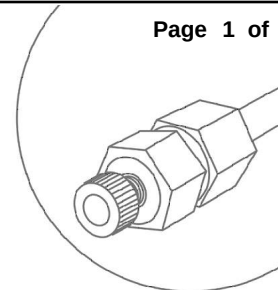

Products used in this Application:

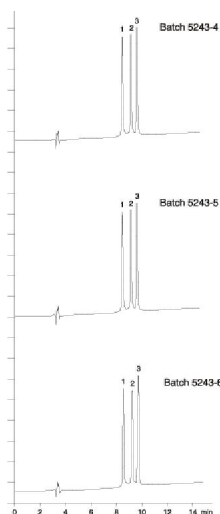

ANALYTES:

- 1 Equine cytochrome c
- 2 Bovine cytochrome c
- 3 Canine cytochrome c

