


Cold Pressed Oragne Oil

Column: Zebron ZB-5ms
Dimensions: 30meter x 0.25mm x 0.25µm
Order No.: 7HG-G010-11
Injection: Split 25:1 @ 250°C, 1µL
Carrier Gas: 1.0mL/min for 3 min to 0.8mL at 0.2mL/min
Oven Program: 40°C for 4 min to 120°C at 8°C/min to 200°C at 15°C/min to 300°C at 25°C/min
Detector: MSD; 35-350amu
Sample:

1. α-Pinene
2. β-Phellandrine
3. β-Myrcene
4. Octanal
5. 3-Carene
6. Limonene
7. Lanilool
8. Nonanal
9. α-Terpineol
10. Decanal

Gold Seal: Phenomenex single groove (AGO-7518)
Liner: Single Taper Gooseneck (AG0-4680)
Septum: Phenomenex PhenoRed-400 (AG0-4690)

Notes: Sample was 10% in methanol

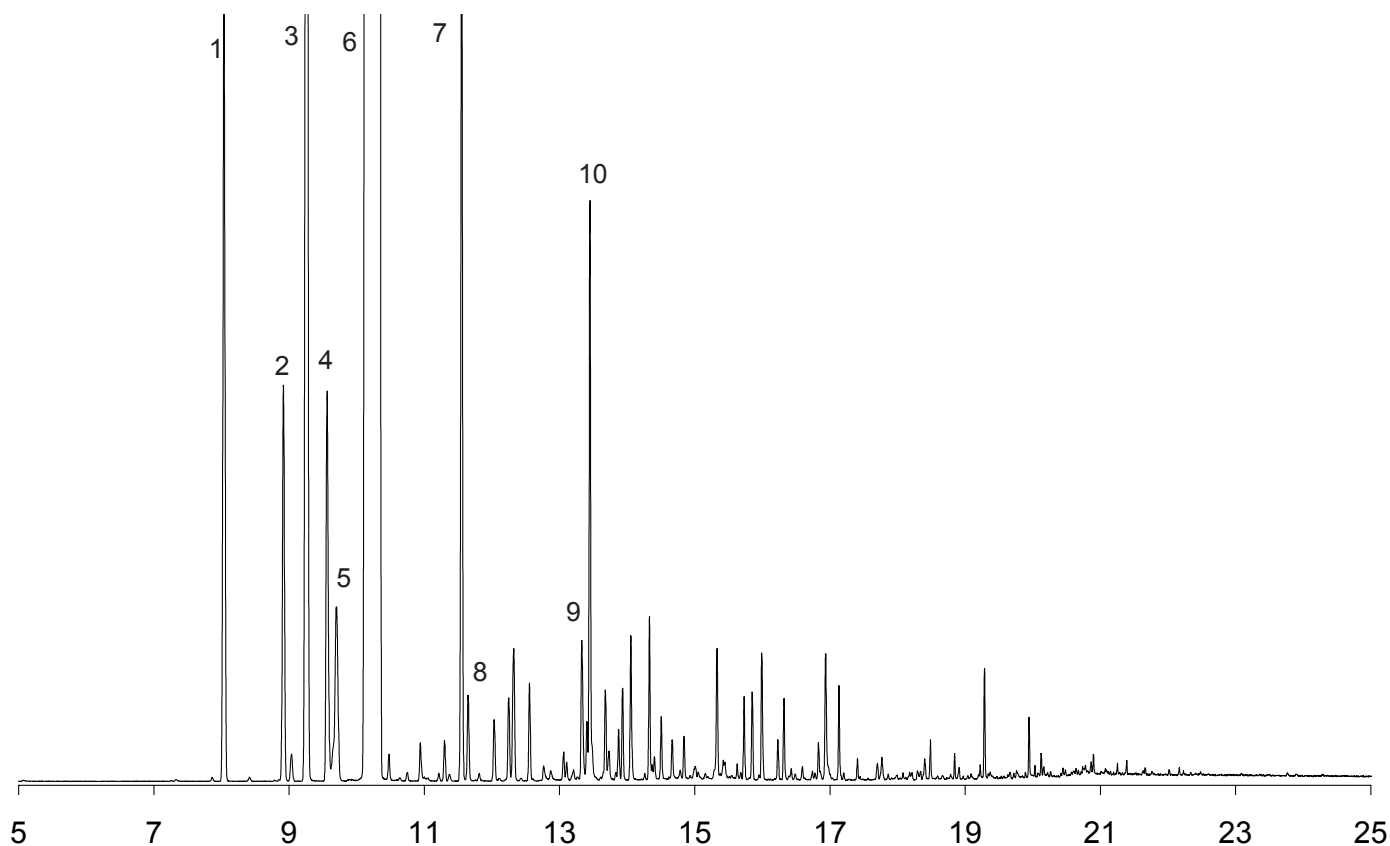

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with traditionSM

USA mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	Puerto Rico 271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926 (800) 541-HPLC (310) 328-7768 info@phenomenex.com	Canada 411 Madrid Ave. Torrance, CA 90501-1430 USA (800) 543-3681 (310) 328-7768 info@phenomenex.com	United Kingdom Queens Avenue, Hurdsfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England 01625-501367 01625-501796 ukinfo@phenomenex.com	Germany Zeppelinstr. 5 63741 Aschaffenburg Deutschland 06021-58830-0 06021-58830-11 anfrage@phenomenex.com	New Zealand P.O. Box 31-601 Milford Auckland New Zealand 09-4780951 09-4780952 info@phenomenex.co.nz	Australia P.O. Box 264 Pennant Hills NSW 1715 Australia 1800-553-929 1800-553-923 info@phenomenex.com.au
--	--	---	---	---	---	--