


Polynuclear Aromatic Hydrocarbons (PAHs)

Column:	Zebtron ZB-5ms	Gold Seal:	Phenomenex single groove (AGO-7518)
Dimensions:	30meter x 0.25mm x 0.25µm	Liner:	Single Taper Gooseneck (AGO-4680)
Order No.:	7HG-G010-11	Septum:	Phenomenex PhenoRed-400 (AGO-4690)
Injection:	Split 15:1 @ 260°C, 1µL		
Carrier Gas:	Constant Flow Helium 1.4 mL/min		
Oven Program:	170°C to 240°C at 15°C/min to 275°C at 4/min to 320°C at 10°C/min for 5 min		
Detector:	MSD; 50-350amu		
Sample:	1. Naphthalene 2. Acenaphthene 3. Acenaphthylene 4. Fluorene 5. Phenanthrene 6. Anthracene 7. Fluoranthene 8. Pyrene 9. Benzo[a]anthracene 10. Chrysene 11. Benzo[b]fluoranthene 12. Benzo[k]fluoranthene 13. Benzo[a]pyrene 14. Indeno[1,2,3-cd]pyrene 15. Dibenz[a,h]anthracene 16. Benzo[g,h,i]perylene		

Notes: All analytes are 100ppm in methylene chloride.

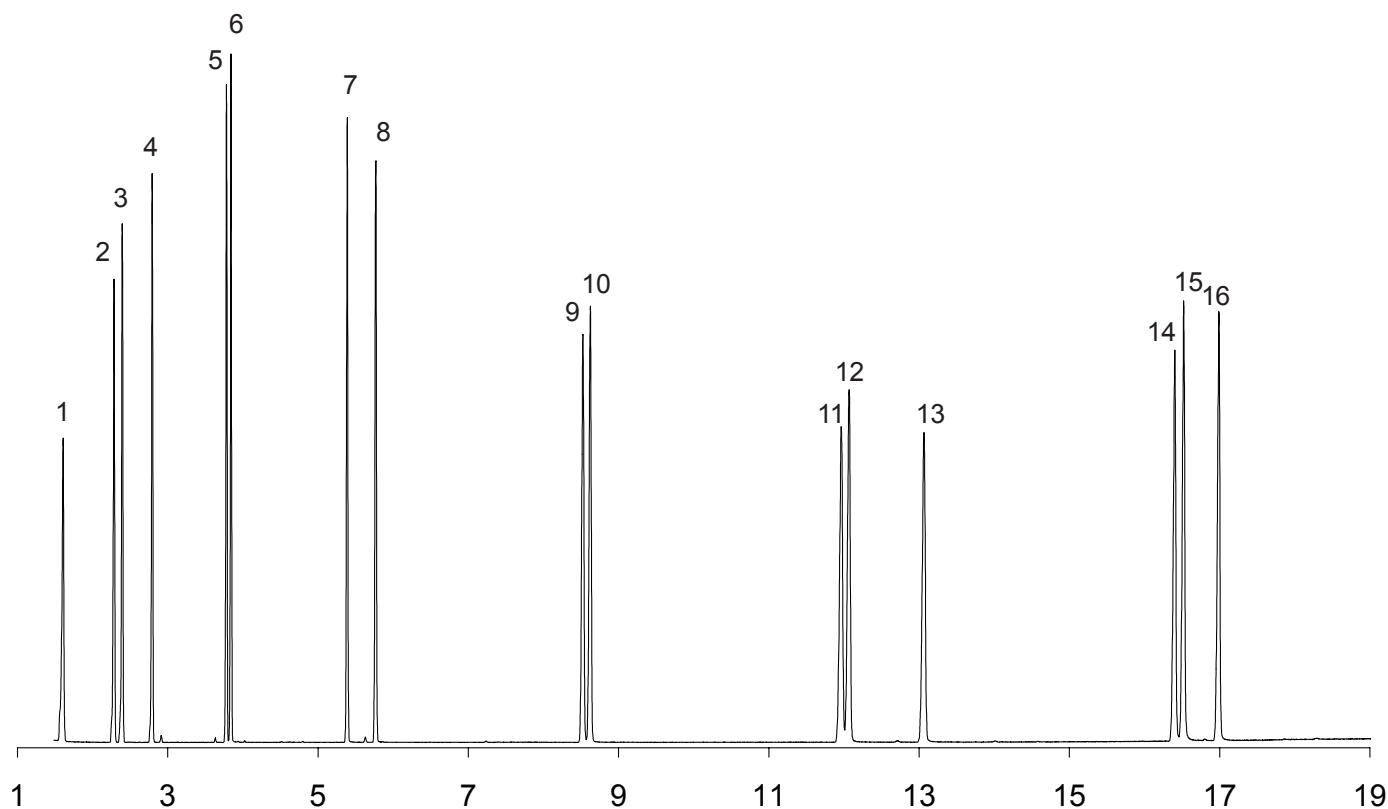

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with traditionSM

USA mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	Puerto Rico 271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926 (800) 541-HPLC (310) 328-7768 info@phenomenex.com	Canada 411 Madrid Ave. Torrance, CA 90501-1430 USA (800) 543-3681 (310) 328-7768 info@phenomenex.com	United Kingdom Queens Avenue, Hurdsfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England 01625-501367 01625-501796 ukinfo@phenomenex.com	Germany Zeppelinstr. 5 63741 Aschaffenburg Deutschland 06021-58830-0 06021-58830-11 anfrage@phenomenex.com	New Zealand P.O. Box 31-601 Milford Auckland New Zealand 09-4780951 09-4780952 info@phenomenex.co.nz	Australia P.O. Box 264 Pennant Hills NSW 1715 Australia 1800-553-929 1800-553-923 info@phenomenex.com.au
--	--	---	---	---	---	--