

RNA Fragments

Column: Jupiter 5u C4 300A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4167-E0

Elution Type: Gradient

Eluent A: 0.1 M K3PO4 in 0.75% isopropanol, pH 6.8


Eluent B: 1.0 M (NH4)SO4 + 0.1M K3PO4 in 0.75% isopropanol, pH 6.8

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	0	100
	2	20	70	30
	3	30	100	0


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 260 nm (ambient)

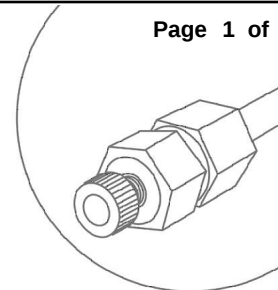

Products used in this Application:

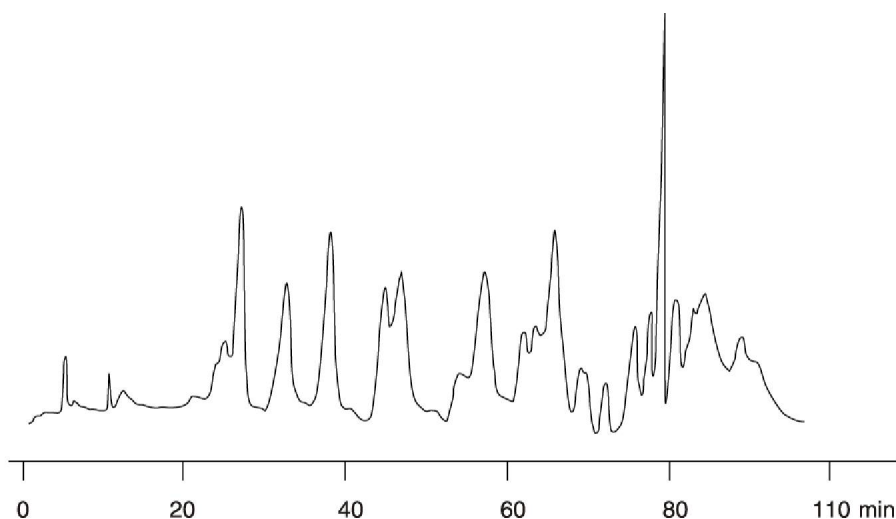Jupiter

ANALYTES:

1 tRNA fragment



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au