


Polynuclear Aromatic Hydrocarbons (PAHs)

Column:	Zebtron ZB-5ms	Inlet Seal:	Phenomenex single groove (AG0-7518)
Dimensions:	30meter x 0.25mm x 0.25µm	Liner:	Single Taper/Goosneck (AG0-4680)
Order No.:	7HG-G010-11	Septum:	Phenomenex PhenoRed-400 (AG0-4690)
Injection:	Split 15:1 @ 285°C, 1µL		
Carrier Gas:	Constant Flow Helium 1.4 mL/min		
Oven Program:	155°C to 250°C at 25°C/min to 280°C at 5°C/min to 320°C at 25°C/min for 2 min		
Detector:	MSD; 40-400 amu		
Sample:	- Naphthalene - Acenaphthene - Acenaphthylene - Fluorene - Phenanthrene - Anthracene - Fluoranthene - Pyrene - Benzo[a]anthracene - Chrysene - Benzo[b]fluoranthene - Benzo[k]fluoranthene - Benzo[a]pyrene - Indeno[1,2,3-cd]pyrene - Dibenz[a,h]anthracene - Benzo[g,h,i]perylene		

Notes: Analytes are at 10ppm in Methylene Chloride.

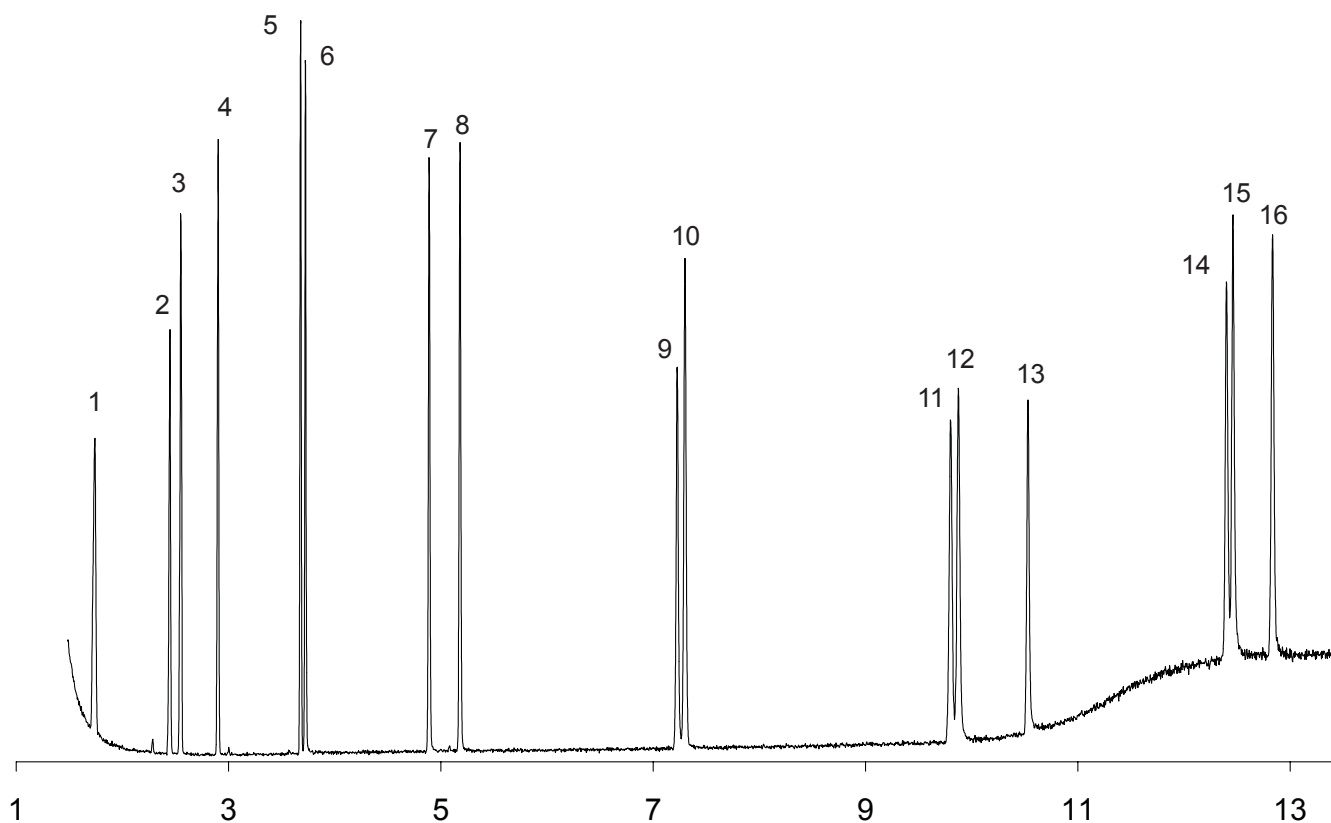

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with traditionSM

USA mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	Puerto Rico 271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926 (800) 541-HPLC (310) 328-7768 info@phenomenex.com	Canada 411 Madrid Ave. Torrance, CA 90501-1430 USA (800) 543-3681 (310) 328-7768 info@phenomenex.com	United Kingdom Queens Avenue, Hurdfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England 01625-501367 01625-501796 ukinfo@phenomenex.com	Germany Zeppelinstr. 5 63741 Aschaffenburg Deutschland 06021-58830-0 06021-58830-11 anfrage@phenomenex.com	New Zealand P.O. Box 31-601 Milford Auckland New Zealand 09-4780951 09-4780952 info@phenomenex.co.nz	Australia P.O. Box 264 Pennant Hills NSW 1715 Australia 1800-553-929 1800-553-923 info@phenomenex.com.au
--	--	---	--	---	---	--