


Benzodiazepines

Column: Zebron ZB-1ms
Dimensions: 12meter x 0.20mm x 0.33µm
Order No.: 7DE-G011-14
Injection: Splitless @ 280°C, 3µL
Carrier Gas: Constant Flow Helium 1.2 mL/min
Oven Program: 150°C to 200°C at 30°C/min to 300 at 15°C/min for 2 min
Detector: MSD; 45-450 amu
Sample:

1. Oxazepam
2. Lorazepam
3. Diazepam
4. Nordazepam
5. Temazepam
6. Flurazepam

Gold Seal: Phenomenex single groove (AG0-7518)
Liner: Split/Splitless Single Taper (AG0-7516)
Septum: Phenomenex PhenoRed-400 (AG0-4690)

Notes: Analytes were 15ng on column

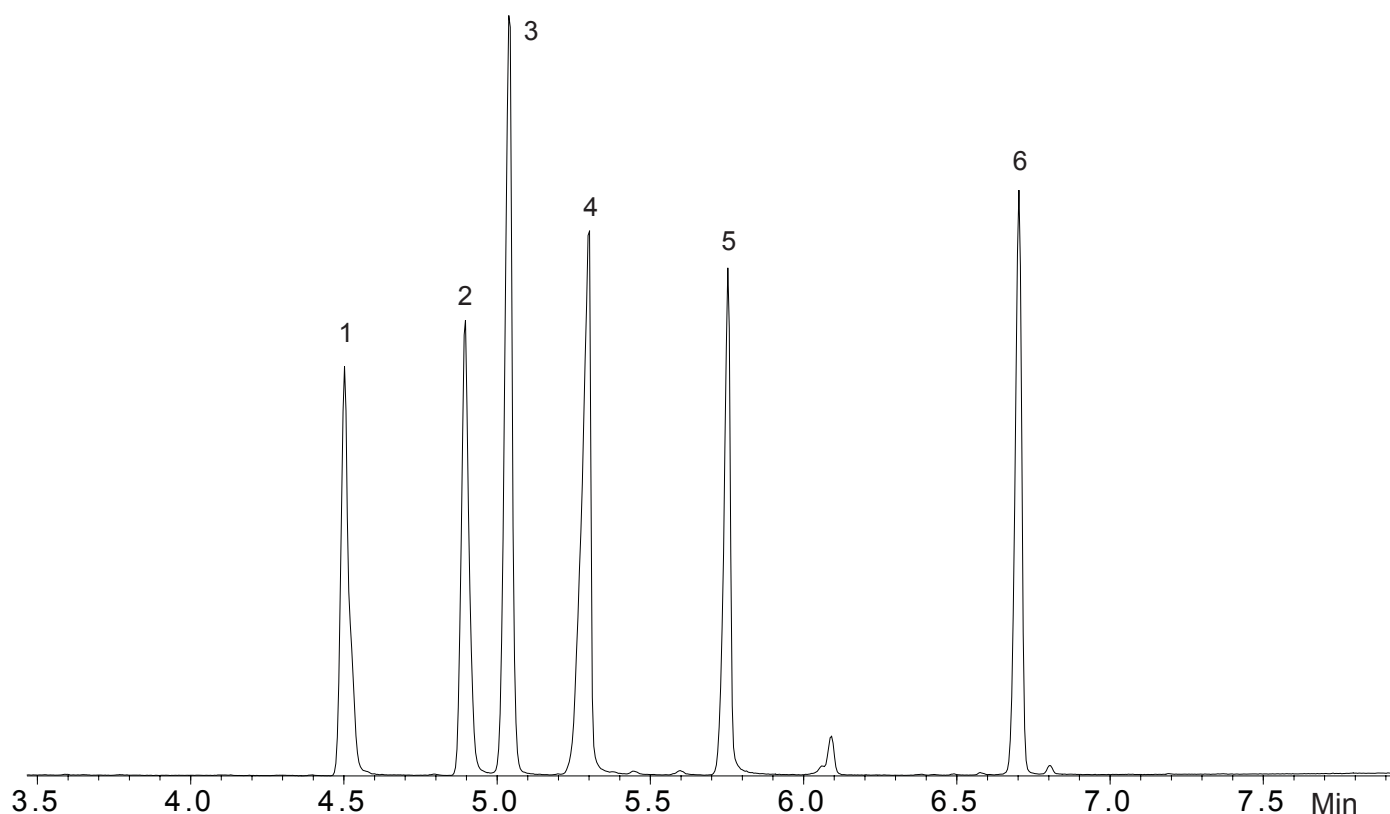

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.


phenomenex
 ...breaking with tradition™

USA
 mail: 411 Madrid Avenue
 Torrance, CA
 90501-1430
 USA
 tel.: 310-212-0555
 fax: 310-328-7768
 email: info@phenomenex.com

Puerto Rico
 273 Sierra Morena,
 Suite #104
 San Juan,
 Puerto Rico 00926
 800-541-HPLC
 310-328-7768
 info@phenomenex.com

Canada
 411 Madrid Avenue
 Torrance, CA
 90501-1430
 USA
 800-543-3681
 310-328-7768
 info@phenomenex.com

France
 118/130
 avenue Jean Jaurès
 75171 Paris Cedex 19
 France
 01 72 72 92 10
 01 72 72 92 11
franceinfo@phenomenex.com

United Kingdom
 Queens Avenue,
 Hurdfield Ind. Est.,
 Macclesfield, Cheshire
 SK10 2BN, England
 01625-501367
 01625-501796
ukinfo@phenomenex.com

Germany
 Zeppelinstr. 5
 63741 Aschaffenburg
 Deutschland
 06021-58830-0
 06021-58830-11
anfrage@phenomenex.com

New Zealand
 P O Box 31601
 Milford
 Auckland
 New Zealand
 09-4780951
 09-4780952
info@phenomenex.co.nz

Australia
 PO Box 4084
 Lane Cove, NSW 2066
 Australia
 02-9428-6444
 02-9428-6445
info@phenomenex.com.au