

Proteins II

Column: Jupiter 5u C4 300A

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4167-E0

Elution Type: Gradient

Eluent A: 0.1% TFA in water


Eluent B: 0.1% TFA in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	1	80	20


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 220 nm (ambient)

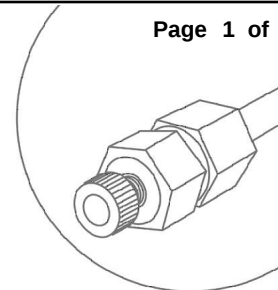

Products used in this Application:

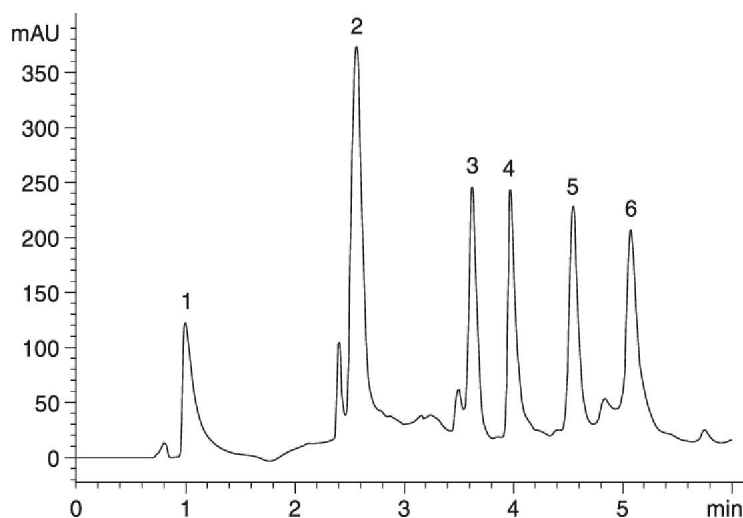

ANALYTES:

- 1 Alkaline phosphatase
- 2 Cyanocobalamin
- 3 RNase
- 4 Insulin
- 5 Transferrin
- 6 Trypsin inhibitor

