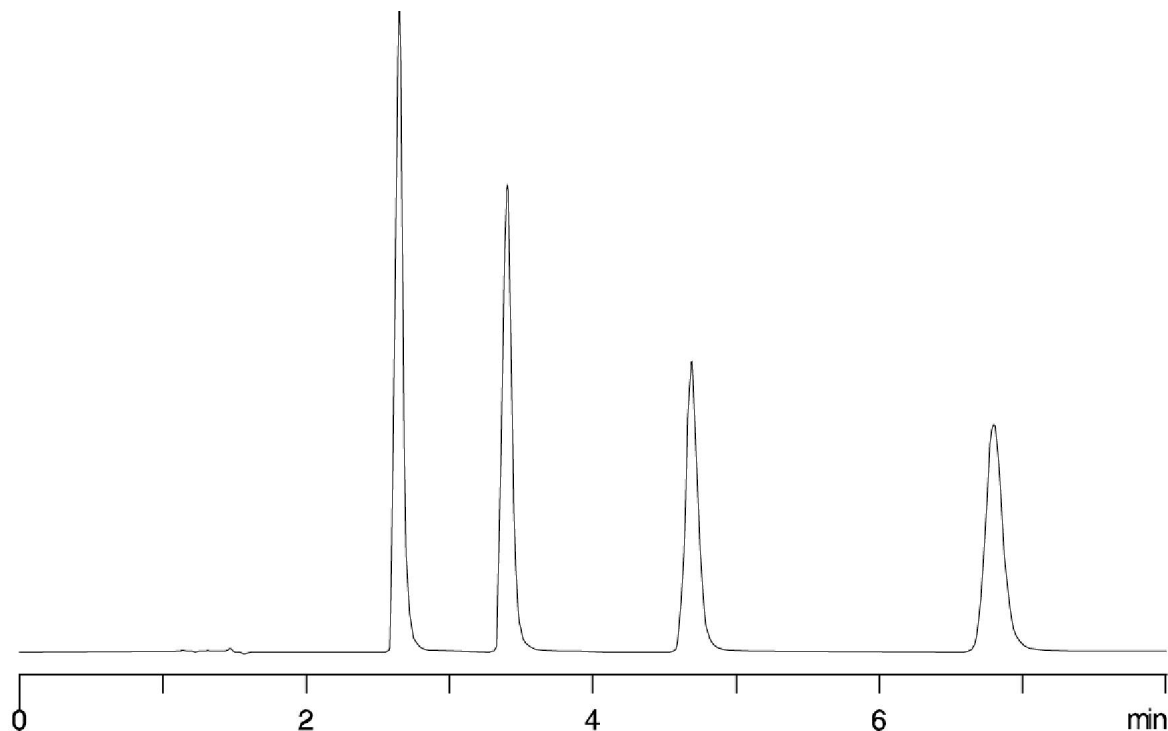


Parabens on Synergi Max-RP

G A L A X Y
S E R I E S**Column:** Synergi 4 μ m C12 Max-RP**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4337-E0**Elution Type:** Isocratic**Eluent A:** Acetonitrile**Pct. A:** 50**Eluent B:** Water**Pct. B:** 50**Flow Rate:** 1 mL/min**Temperature:** 30 °C**Detection:** UV-Vis Abs. - Variable Wave. @ 254 nm**SAMPLE ANALYTES:**

- 1 Methyl paraben
- 2 Ethyl paraben
- 3 Propyl paraben
- 4 Butyl paraben

SPE:**For more information call your Phenomenex Representative at (310) 212-0555**

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 12642/BRO

Page 1 of 1