


APPLICATION NOTE - HPLC
Parabens on C18

LUNA™

G A L A X Y™
S E R I E S

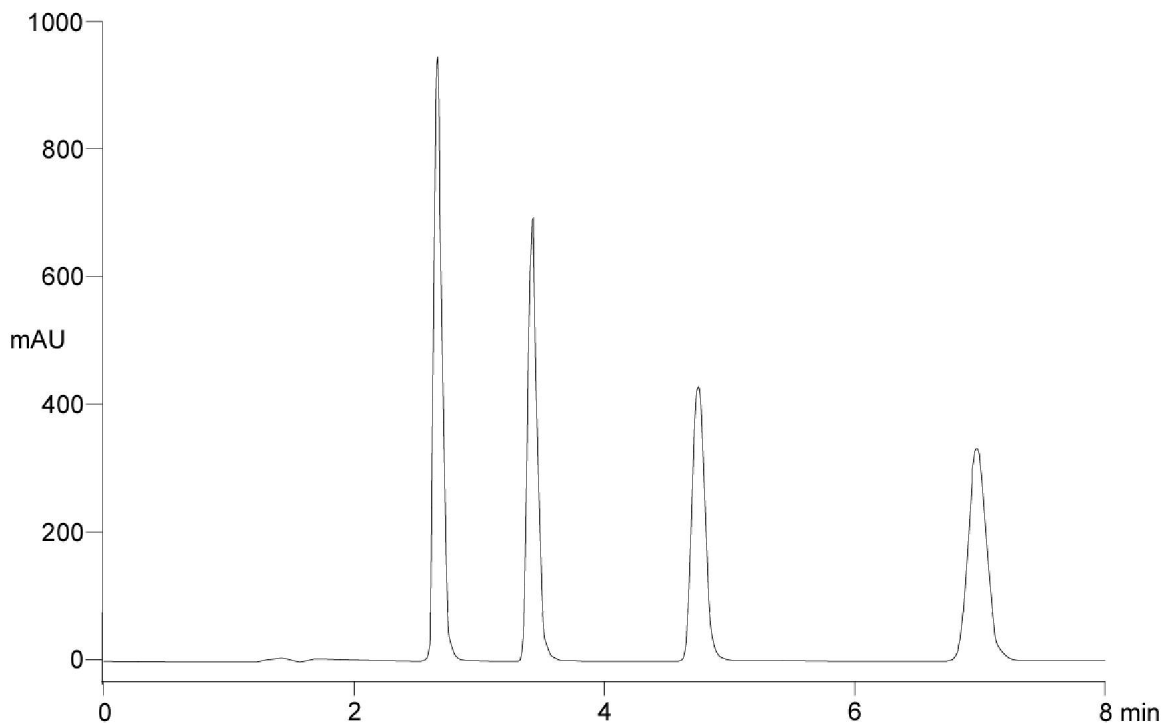

Column: Luna 5 µm C18 / ODS (2)

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Pct. A: 50

Eluent B: Water

Pct. B: 50

Flow Rate: 1 mL/min

Temperature: 30 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTES:

- 1 Methyl paraben
- 2 Ethyl paraben
- 3 Propyl paraben
- 4 Butyl paraben

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 12643

Page 1 of 1