

Echinacea (Echinacea angustifolia)

Column: Luna 5u C18(2) 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Isocratic

Eluent A: Acetonitrile/water, 0.1% phosphoric Acid (10:90)

Pct. A: 50

Eluent B: Acetonitrile/water, 0.1% phosphoric Acid (25:75)


Pct. B: 50

Flow Rate: 1.5 mL/min

Col. Temp.: ambient


Detection: Unknown (0 °C)

Page 1 of 1

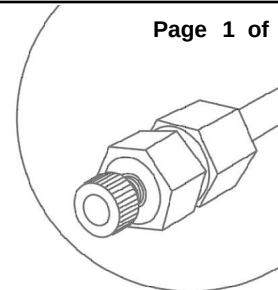

Products used in this Application:

LUNA®

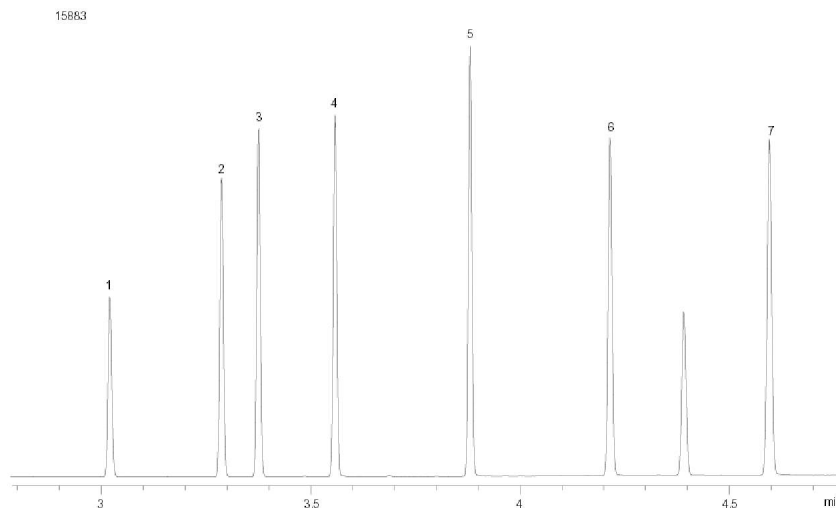

ANALYTES:

- 1 Caffeic acid
- 2 Chlorogenic acid
- 3 Echinacoside
- 4 Cichoric acid

