

Dextran

Column: PolySep-GFC-P 4000

Dimensions: 300 x 7.8 mm ID

Order No: 00H-3144-K0

Elution Type: Isocratic

Eluent A: Water



Pct. A: 100

Flow Rate: 0.8 mL/min

Col. Temp.: ambient

Detection: Refractive Index (RI) (ambient)

Page 1 of 1

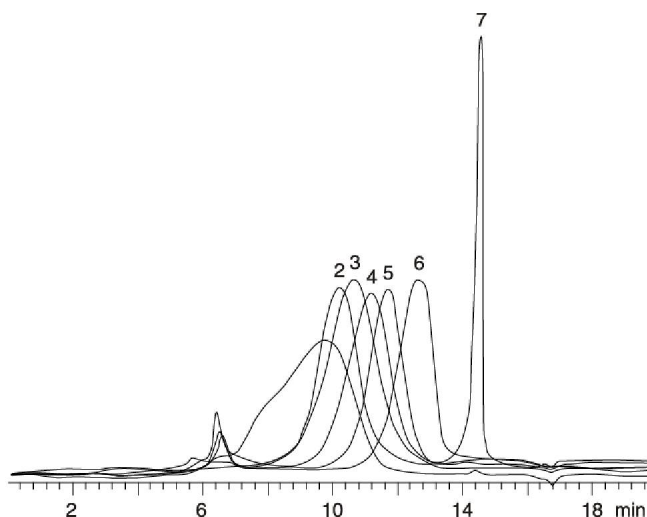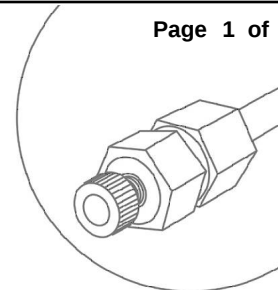

ANALYTES:

- 1 Dextran 500 KD
- 2 Dextran 110 KD
- 3 Dextran 70 KD
- 4 Dextran 40 KD
- 5 Dextran 15 KD
- 6 Dextran 6 KD
- 7 Ethylene glycol 62.07 D

