


Aflatoxins Chromatographic Comparison

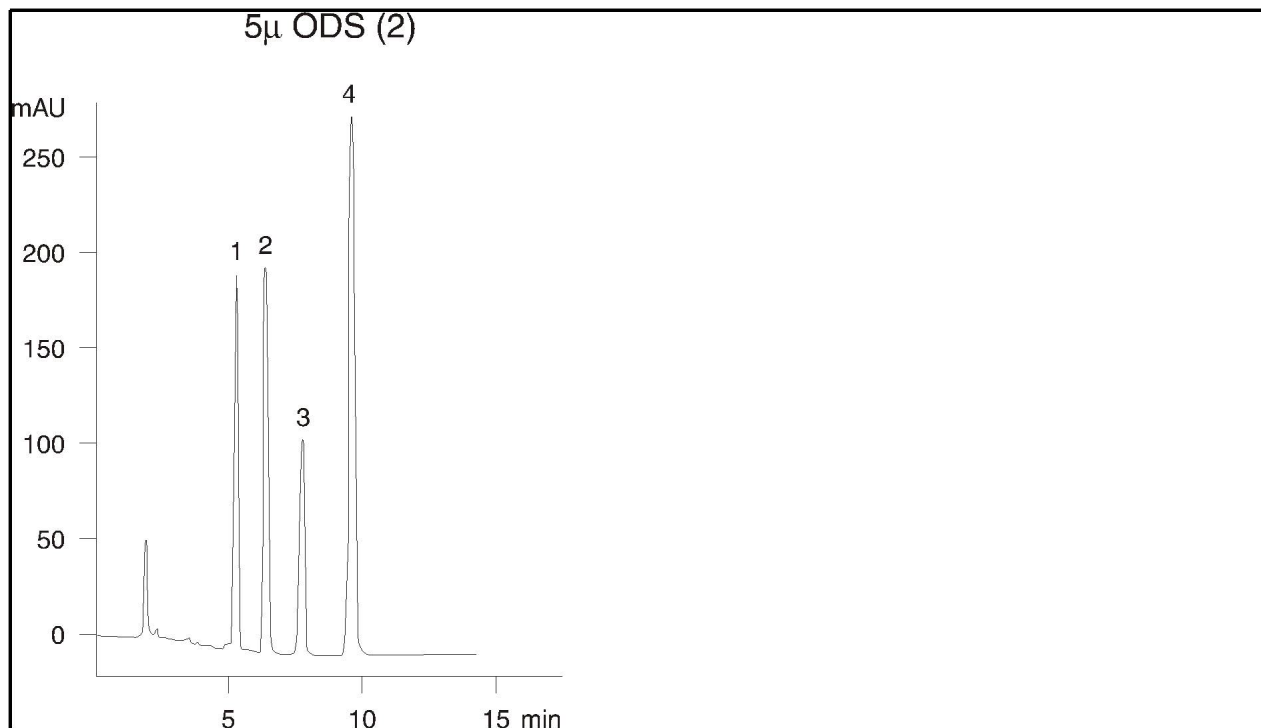

Column: Prodigy 5 μ m C18 / ODS (2)

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Water

Pct. A: 55

Eluent B: Methanol

Pct. B: 45

Flow Rate: 0.8 mL/min

Temperature: 40 °C

Detection: UV-Vis Abs. - Variable Wave. @ 365 nm

SAMPLE ANALYTE

- 1 Aflatoxin G2
- 2 Aflatoxin G1
- 3 Aflatoxin B2
- 4 Aflatoxin B1

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 5478/CAT99/227