



Amino Acids - Underivatized

Column: Aqua 5u C18 125A**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4299-E0**Elution Type:** Isocratic**Eluent A:** Water, containing 0.3% heptafluorobutyric acid**Pct. A:** 75**Eluent B:** Methanol, containing 0.3% heptafluorobutyric acid**Pct. B:** 25**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) (ambient)**Comment:** Positive APCI-MS (Finnigan LCQ)

Page 1 of 1

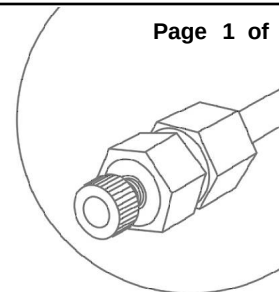

Products used in this Application:

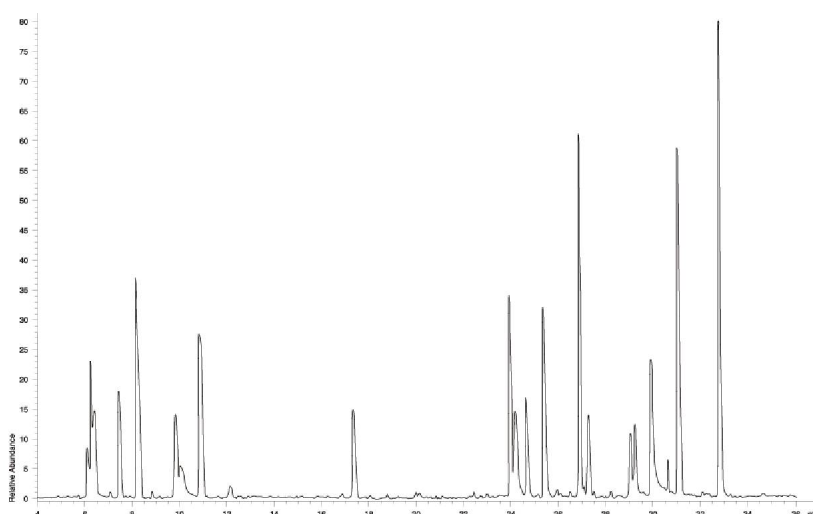AQUA™

ANALYTES:

- | | |
|-----------------|------------------|
| 1 Glycine | 17 Isoleucine |
| 2 Asparagine | 18 Leucine |
| 3 Serine | 19 Phenylalanine |
| 4 Aspartic acid | 20 Tryptophan |
| 5 Glutamine | |
| 6 Threonine | |
| 7 Alanine | |
| 8 Glutamic acid | |
| 9 Cysteine | |
| 10 Proline | |
| 11 Lysine | |
| 12 Valine | |
| 13 Histidine | |
| 14 Methionine | |
| 15 Tyrosine | |
| 16 Arginine | |

©Copyright 2000-2006 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.phenomenex®
...breaking with tradition™**USA**
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com**Puerto Rico**
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com**Canada**
(800) 543-3681
(310) 328-7768
info@phenomenex.com**France**
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com**United Kingdom**
01625-501367
01625-501796
ukinfo@phenomenex.com**Germany**
06021-58830-0
06021-58830-11
anfrage@phenomenex.com**New Zealand**
09-4780951
09-4780952
info@phenomenex.co.nz**Australia**
02-9428-6444
02-9428-6445
info@phenomenex.com.au