



EZ:faast Amino Acid Analysis of Grape Juice (GC/MS)

Page 1 of 1

Column: EZ:faast GC
Phase: Zebron ZB-AAA (proprietary)
Dimensions: 10 meters x 0.25 mm x 0 µm
Order No: KG0-7166
Oven Profile: 30°C/min from 110° to 320°C, hold at 320°C for 1 minute
Carrier Gas: Constant Flow Helium, 1.1 mL/min
Injection: Split 15:1 2.5 µL @ 250 °C
Detection: Mass Spectrometer (MS) (300 °C)

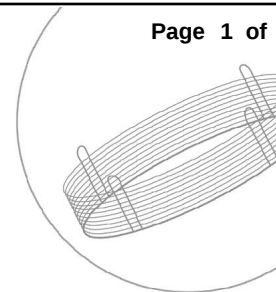

Products used in this Application:

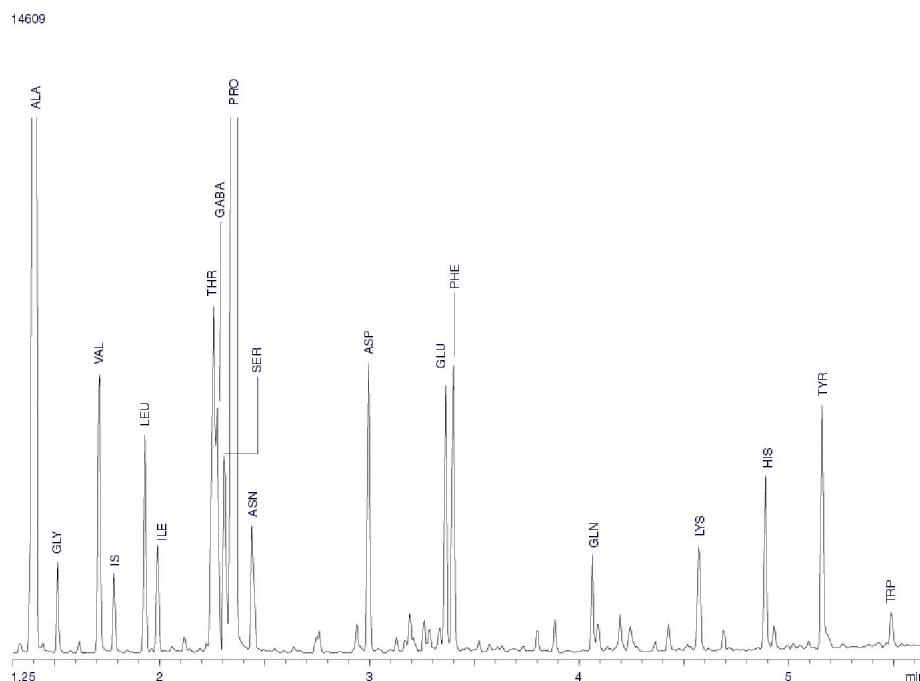EZ:faast[™]**Analyst Note:** Hot needle injection**ANALYTES:**

1 ALA	10 PRO	19 TRP
2 GLY	11 ASN	
3 VAL	12 ASP	
4 IS	13 GLU	
5 LEU	14 PHE	
6 ILE	15 GLN	
7 THR	16 LYS	
8 GABA	17 HIS	
9 SER	18 TYR	

©Copyright 2000-2006 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.phenomenex[®]
...breaking with tradition[™]

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au