

Honey

Column: Rezex RCM-Monosaccharide Ca+ (8%)

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0130-K0

Elution Type: Isocratic

Eluent A: Water


Flow Rate: 0.6 mL/min

Col. Temp.: 75 °C

Detection: Refractive Index (RI) (ambient)


Pct. A: 100

Page 1 of 1

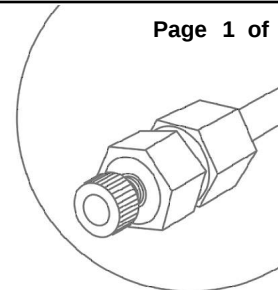

Products used in this Application:

REZEX™

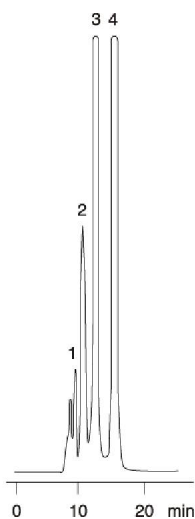

ANALYTES:

- 1 dp 3
- 2 dp 2
- 3 Glucose
- 4 Fructose



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.