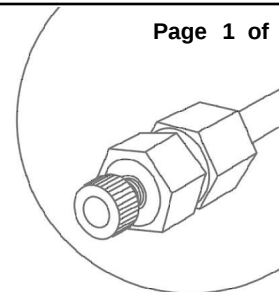



Oligosaccharide Std. (AL0-3038)

Column: Rezex RSO-Oligosaccharide Ag+ 4%**Dimensions:** 200 x 10 mm ID**Order No:** 00P-0133-N0**Elution Type:** Isocratic**Eluent A:** Water**Pct. A:** 100**Flow Rate:** 0.3 mL/min**Col. Temp.:** 85 °C**Detection:** Refractive Index (RI) (85 °C)

Page 1 of 1

**Products used in this Application:****REZEX[™]**

Analyst Note: This standard is used to test the following HPLC columns: REZEX RSO (Ag Oligosaccharide) and RNO (Na Oligosaccharide)
NOTE: Separation conditions and results may vary for your column. Contact your Technical Representative.

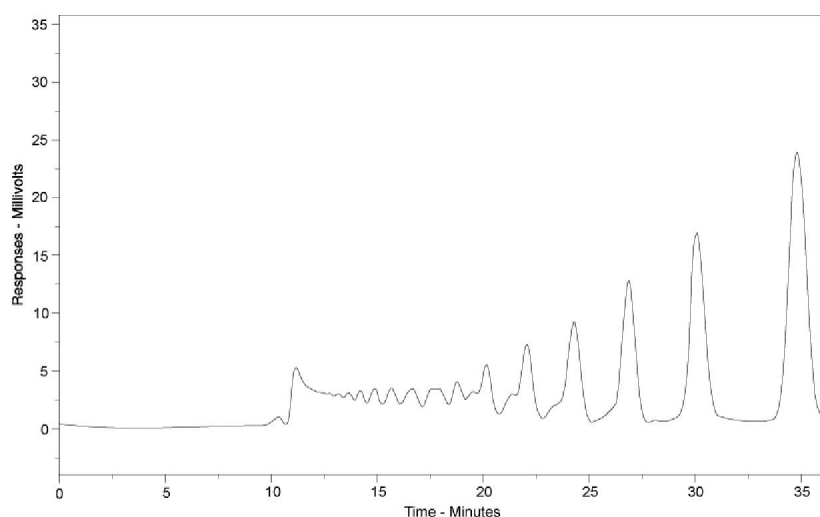

ANALYTES:

- 1 Light corn syrup

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.