

Fluoxetine at pH 2.5 using Onyx C18 3.0mm ID

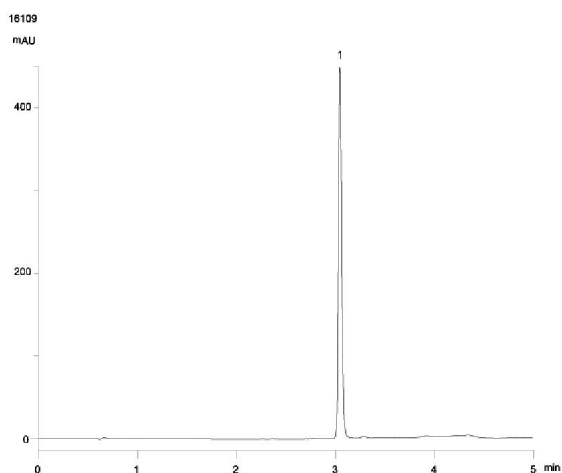
Page 1 of 1

Column: Onyx Monolithic C18**Dimensions:** 100 x 3 mm ID**Order No:** CH0-8158**Elution Type:** Gradient**Eluent A:** 20mM Potassium phosphate pH 2.5 in water**Eluent B:** 20mM Potassium phosphate pH 2.5 in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	3	3	20	80

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 230 nm (ambient)

Products used in this Application:

ONYX[™]

ANALYTES:

1 Fluoxetine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.