


Echinacea (Echinacea purpurea)**Column:** Luna 5u C18(2) 100A**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4252-E0**Elution Type:** Gradient**Eluent A:** Water/acetonitrile (90.5:9.5) + 1% Acetic acid**Eluent B:** Acetonitrile


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	30	60	40
	3	40	100	0
	4	50	100	0

Flow Rate: 1 mL/min**Col. Temp.:** 25 °C**Detection:** UV-Vis Abs. - Variable Wave. @ 284 nm (ambient)

Page 1 of 1

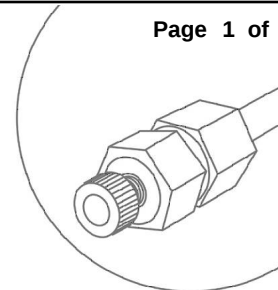

Products used in this Application:

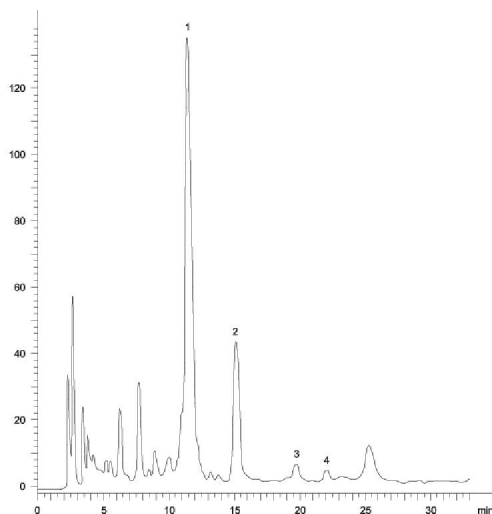LUNA[®]**ANALYTES:**

- 1 Caftaric acid
- 2 Chlorogenic acid
- 3 Caffeic acid
- 4 Syringic acid



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.