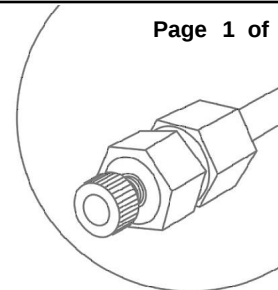
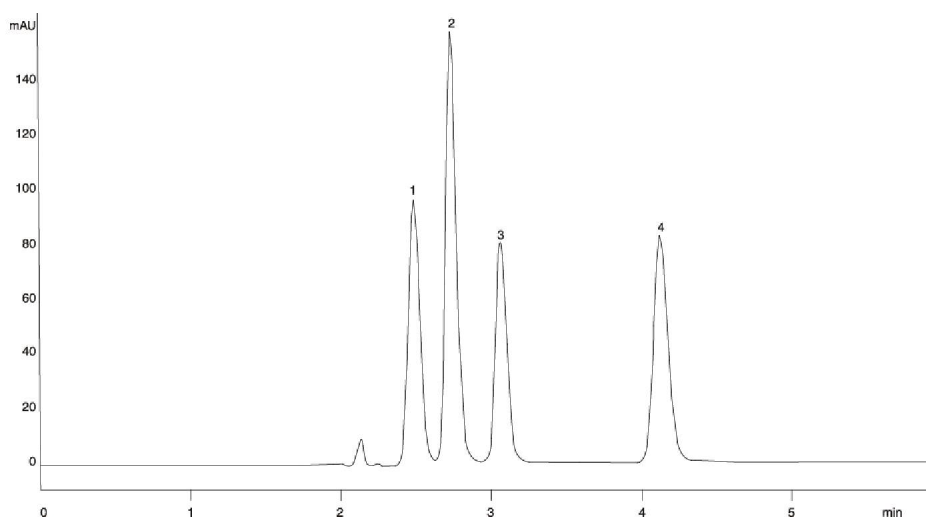


Polyphenols

Column: Luna 5u CN 100A**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4255-E0**Elution Type:** Isocratic**Eluent A:** 5mM Acetic acid**Pct. A:** 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 280 nm (ambient)

Page 1 of 1

**Products used in this Application:****LUNA®****ANALYTES:**

- 1 Pyrogallol
- 2 Hydroquinone
- 3 Resorcinol
- 4 Phenol

**www.phenomenex.com**

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.