

Catecholamines metabolites using Synergi Fusion-RP

Page 1 of 1

Column: Synergi 4u Fusion-RP 80A
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4424-E0
Elution Type: Gradient
Eluent A: 20mM Ammonium formate, pH 3.5
Eluent B: Acetonitrile


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	10	50	50

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 260 nm (ambient)

Products used in this Application:

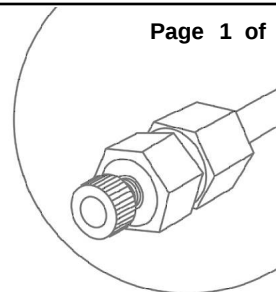synergi[™]

15739

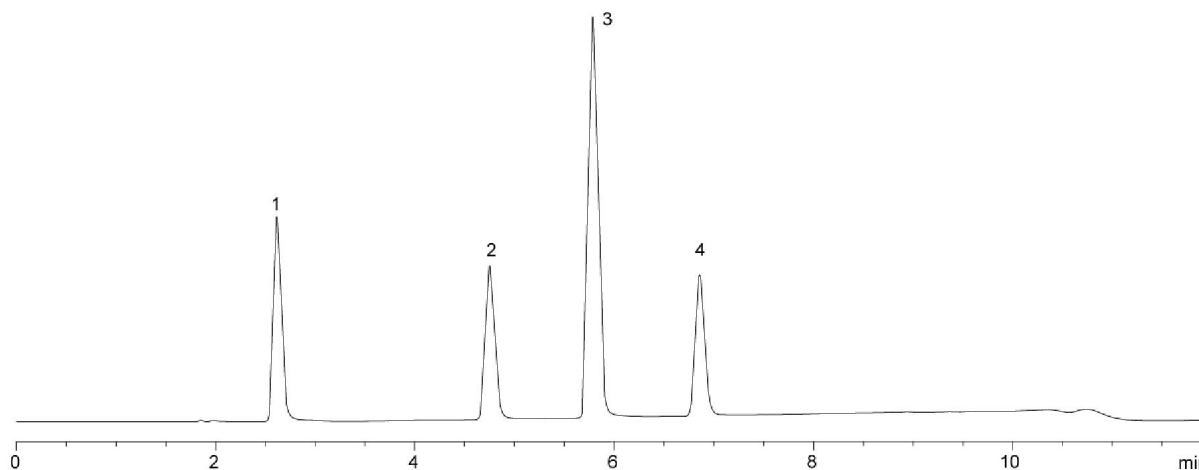

ANALYTES:

- 1 Vanillylmandelic acid
- 2 3,4-Dihydroxyphenylacetic acid (DOPAC)
- 3 5-Hydroxyindoleacetic acid
- 4 Homovanillic acid



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.