


Polar Aromatic Organic Acids by Gemini 5u C6-Phenyl

Page 1 of 1

Column: Gemini 5u C6-Phenyl 110A
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4444-E0
Elution Type: Isocratic
Eluent A: 20mM Potassium phosphate buffer, pH 2.5
Eluent B: Methanol
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Pct. A: 55

Pct. B: 45

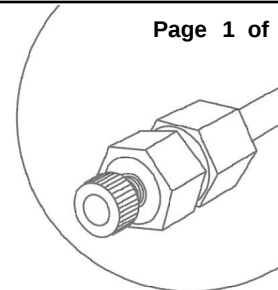


Products used in this Application:



Analyst Note: SYMM (salicylic acid) = 0.66
 Backpressure = 150 bar

15993

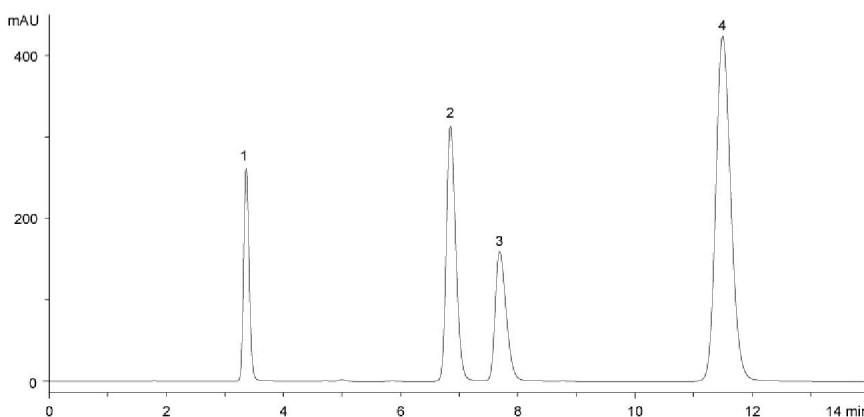

ANALYTES:

- 1 p-Hydroxybenzoic acid
- 2 Benzoic acid
- 3 Salicylic acid
- 4 p-Toluic acid



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.