

Peptides on Synergi Hydro-RP

Page 1 of 1

Column: Synergi 4u Hydro-RP 80A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4375-E0

Elution Type: Gradient

Eluent A: Water with 0.1% TFA

Eluent B: Acetonitrile with 0.1% TFA

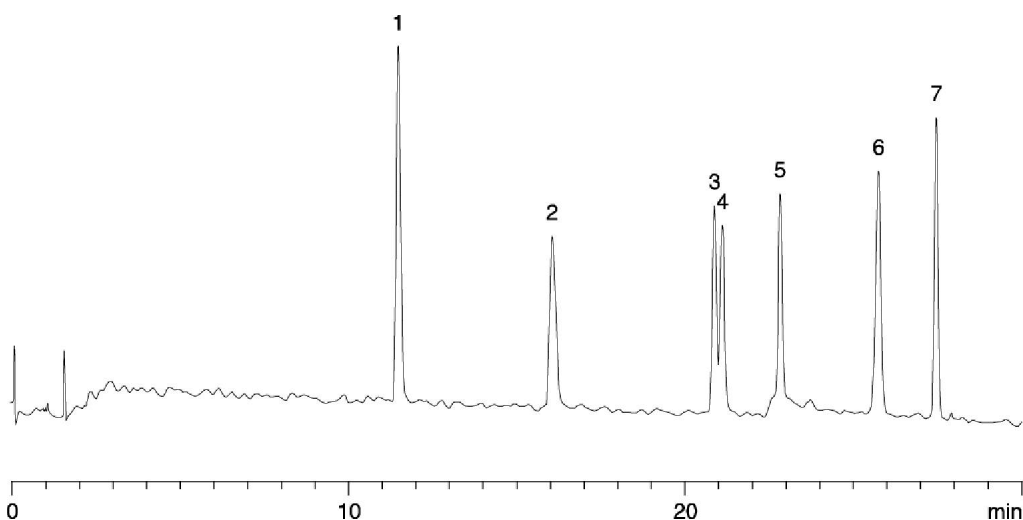
Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	30	55	45

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Products used in this Application:

synergi[™]

ANALYTES:

- 1 Asp-Phe methyl ester
- 2 Bradykinin
- 3 Human Angiotensin I
- 4 Frog Angiotensin I
- 5 Substance P
- 6 Insulin
- 7 N-Formyl-Met-Leu-Phe



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.