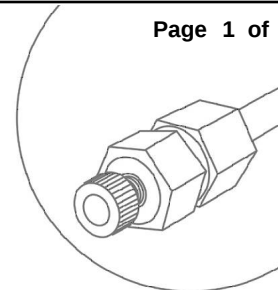
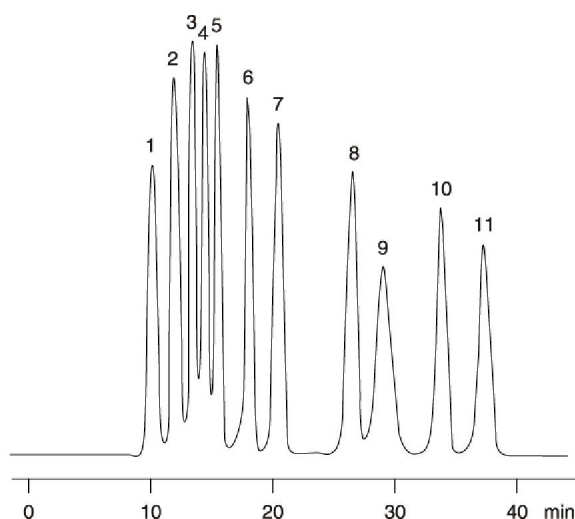


Saccharides

Column: Rezex RPM-Monosaccharide Pb++(8%)**Dimensions:** 300 x 7.8 mm ID**Order No:** 00H-0135-K0**Elution Type:** Isocratic**Eluent A:** Water**Pct. A:** 100**Flow Rate:** 0.6 mL/min**Col. Temp.:** 75 °C**Detection:** Refractive Index (RI) (ambient)

Page 1 of 1

**Products used in this Application:****REZEX™****ANALYTES:**

- 1 Stachyose
- 2 Maltose
- 3 Glucose
- 4 Xylose
- 5 Galactose
- 6 Fructose
- 7 Meso-erythritol
- 8 Mannitol
- 9 Salicin
- 10 Xylitol
- 11 Sorbitol

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.