


Steroids

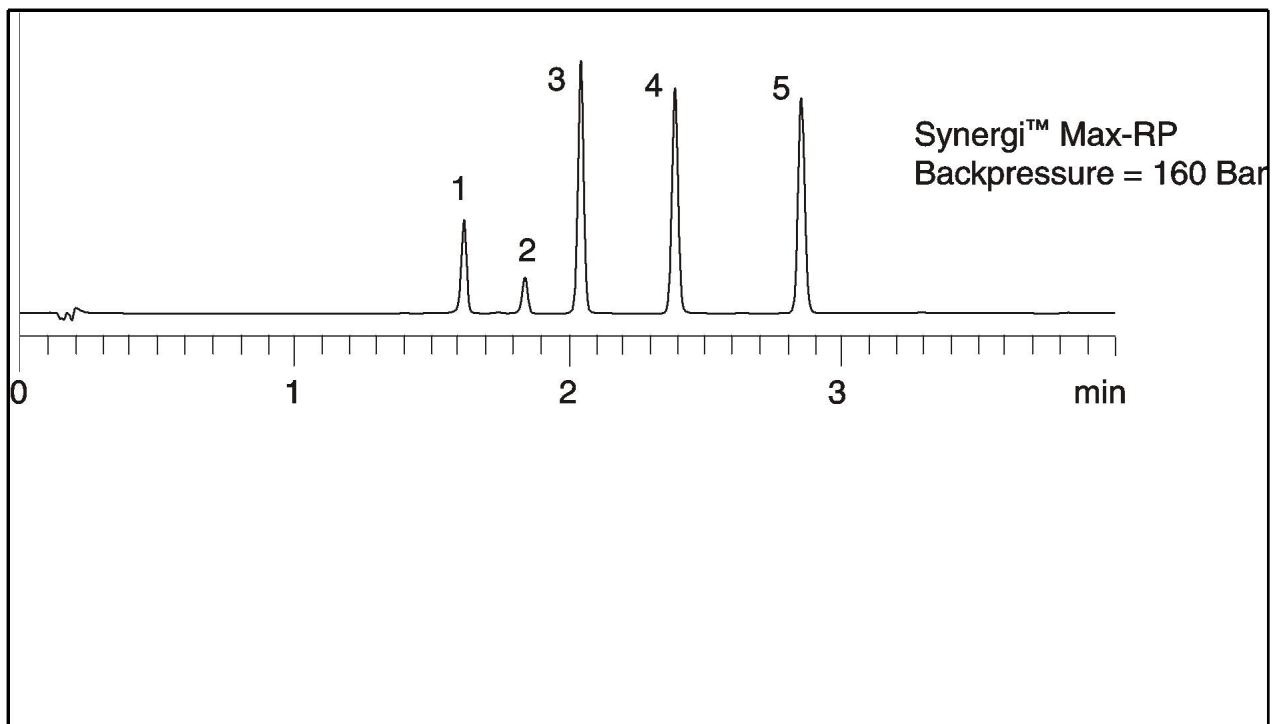

Column: Synergi 4 μ m C12 Max-RP
Dimensions: 50 x 4.6 mm ID
Order No:
Elution Type: Gradient
Eluent A: Water with 0.1% Formic acid
Eluent B: Acetonitrile with 0.1% Formic acid

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	95	5
2	5	5	95

Flow Rate: 4 mL/min
Temperature: 30 °C
Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SPE:**SAMPLE ANALYTE**

- 1 Dexamethasone
- 2 Hydrocortisone
- 3 11-a-Hydroxyprogesterone
- 4 17-a-Hydroxyprogesterone
- 5 Progesterone

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 12480/BRO