


Malic acid

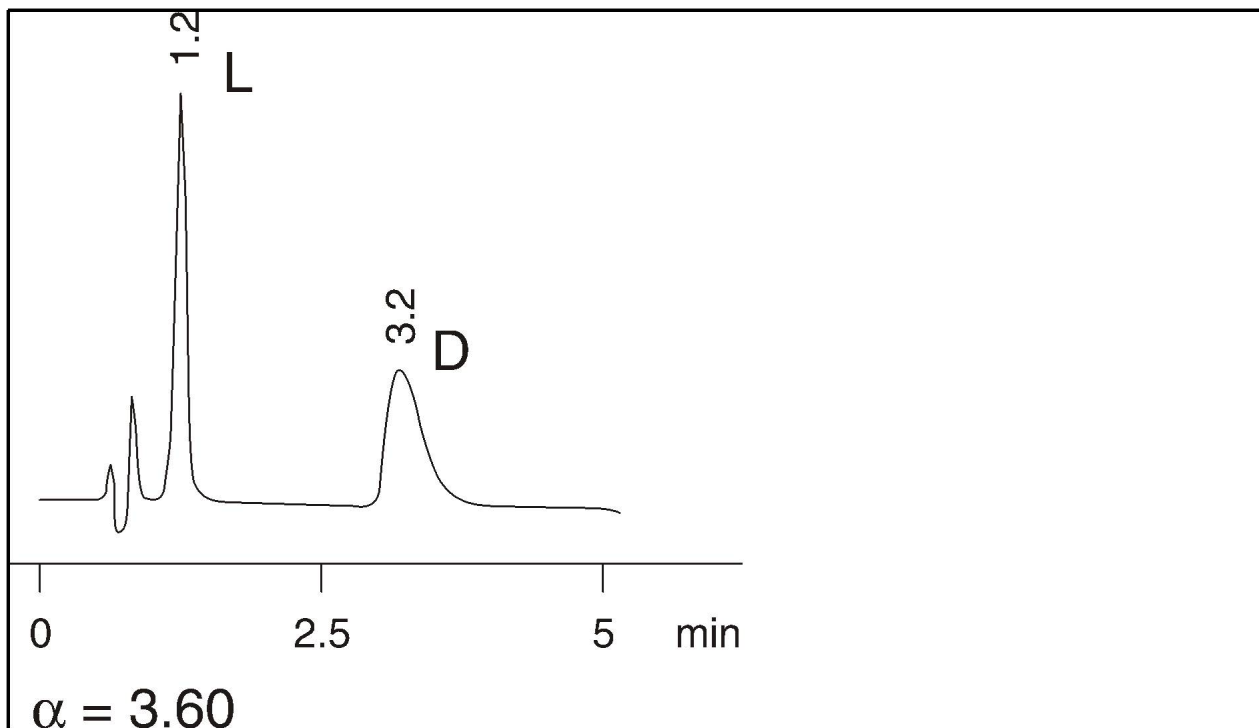

Column: Chirex 5 μ m Chiral IV (ligand exchange) 3126 (R)-PA / D-PA

Dimensions: 50 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 1 mM Copper (II) acetate and 0.1 M ammonium acetate in water **Pct. A:** 85

Eluent B: Isopropanol [pH 4.5] **Pct. B:** 15

Flow Rate: 1 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTE

- 1 L-Malic acid
- 2 D-Malic acid

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14099/CHX/H-05