


Fatty Acids

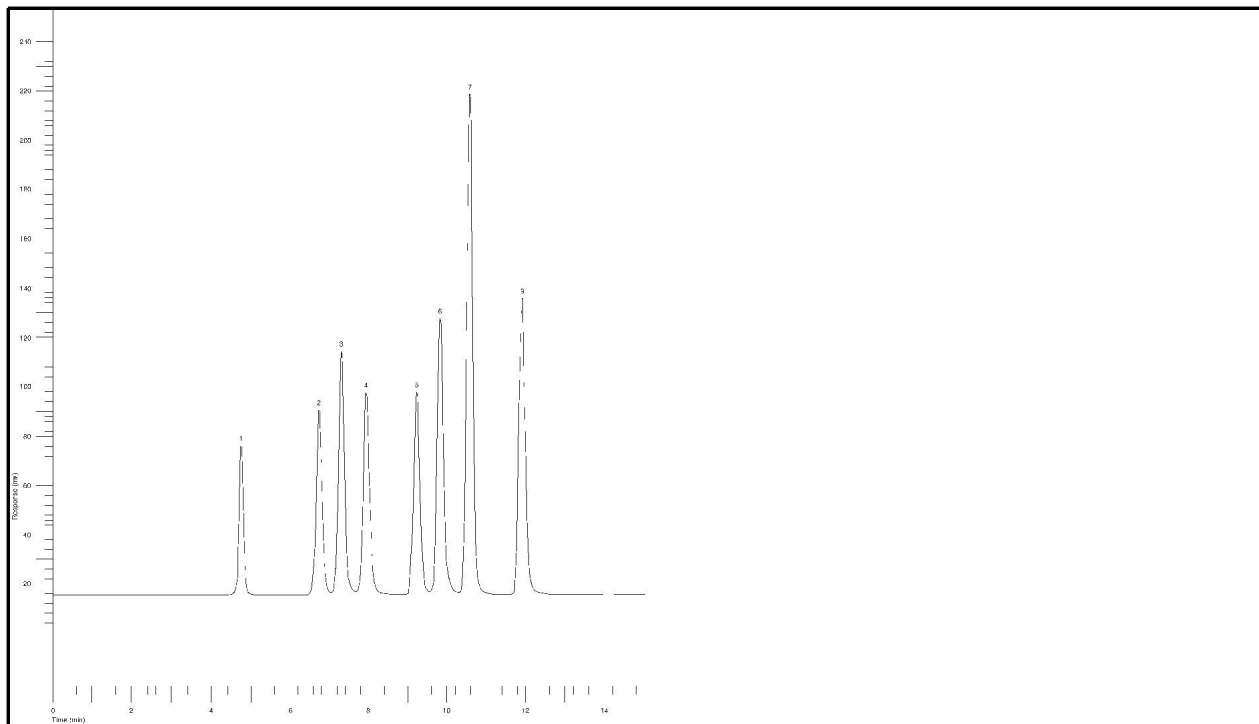

Column: Luna 5 μ m C8 (2)
Dimensions: 150 x 4.6 mm ID
Order No:
Elution Type: Gradient
Eluent A: Acetonitrile
Eluent B: Water (18Mohms DI)

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	70	30
2	10	90	10
3	12	70	30
4	16	70	30

Flow Rate: 0.3 mL/min
Temperature: 22 °C
Detection: Evaporative Light Scattering @

SPE:

SAMPLE ANALYTE

- | | |
|--------------------|----------------------|
| 1 Lauric acid | 7 Heptadecanoic acid |
| 2 Myristic acid | 8 Stearic acid |
| 3 Palmitoleic acid | |
| 4 Linoleic acid | |
| 5 Palmitic acid | |
| 6 Oleic acid | |

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1258/ALG/120