

Sugar Alcohols

Column: Rezex RPM-Monosaccharide Pb++(8%)

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0135-K0

Elution Type: Isocratic

Eluent A: Water


Pct. A: 100

Flow Rate: 0.6 mL/min

Col. Temp.: 75 °C


Detection: Refractive Index (RI) (ambient)

Page 1 of 1

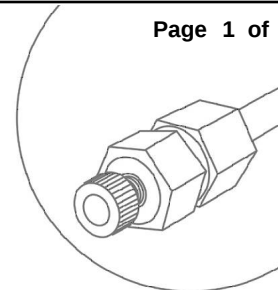

Products used in this Application:

REZEX[™]

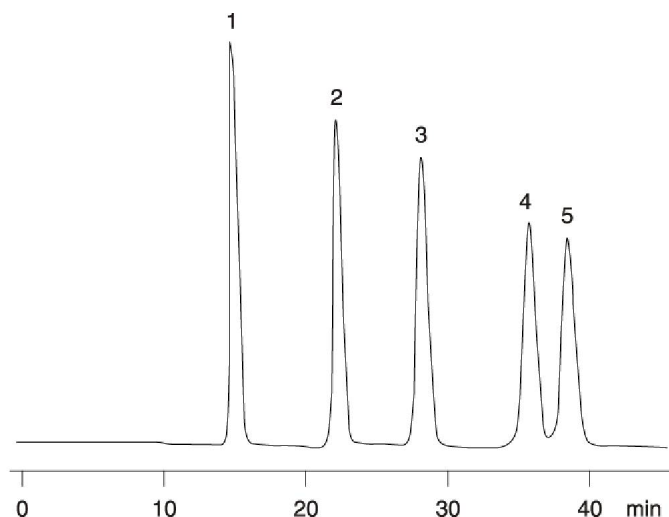

ANALYTES:

- 1 Glucose
- 2 Erythritol
- 3 Mannitol
- 4 Dulcitol
- 5 Sorbitol



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.