


Serine

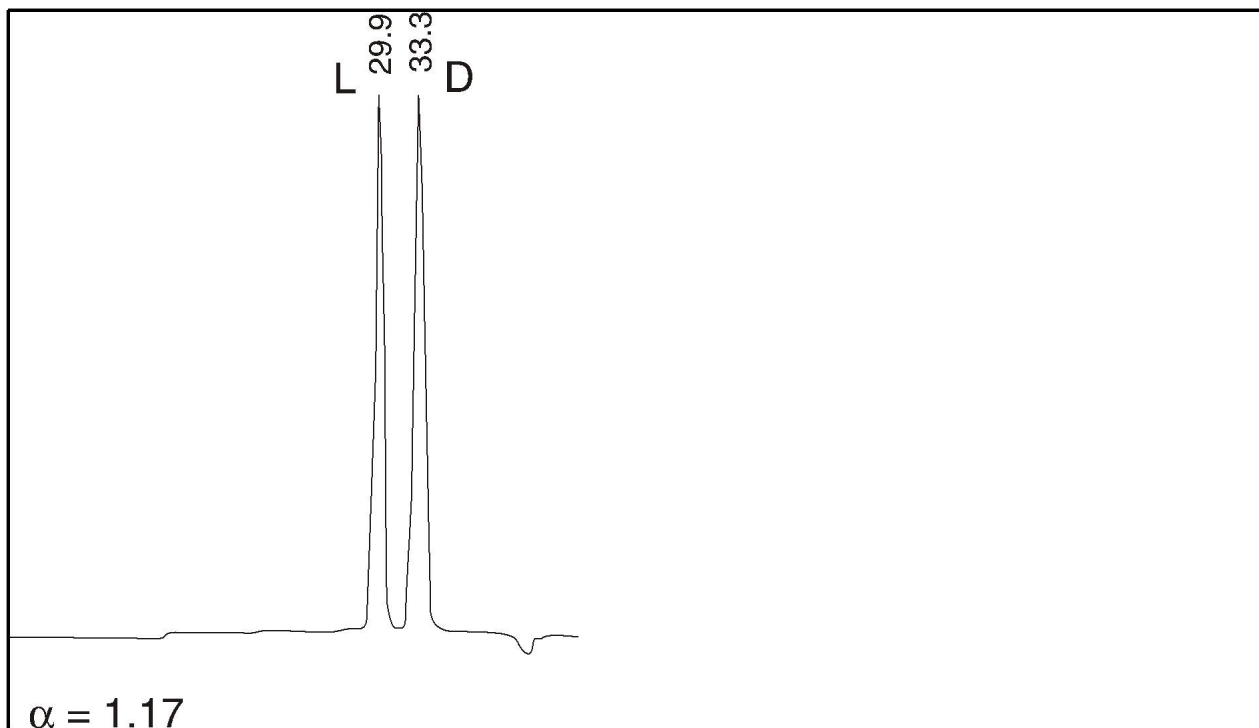

Column: Chirex 5 μ m Chiral IV (ligand exchange) 3126 (R)-PA / D-PA

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 1 mM Copper (II) sulfate in water

Pct. A: 100

Flow Rate: 0.2 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTE

1 L-Serine

2 D-Serine

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14016/CHX/G-05