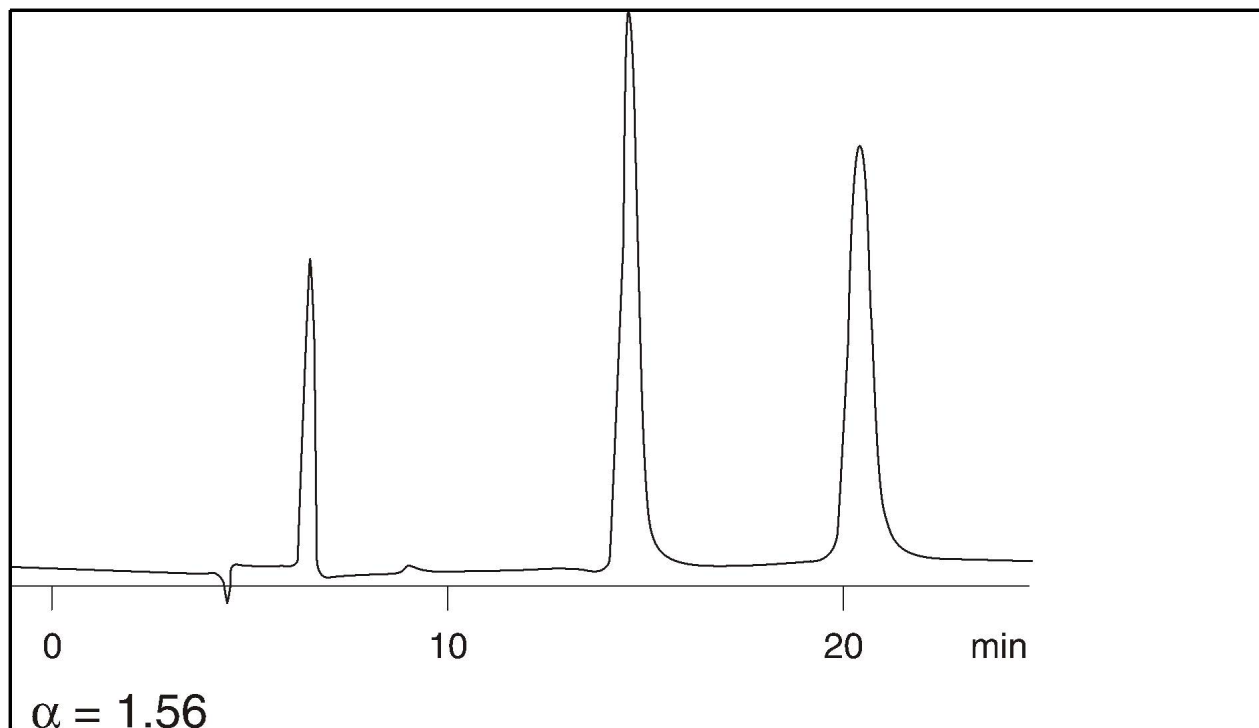


1,1'-Bi-2-naphthol

Column: Chirex 5 μ m Chiral I (Pirkle) 3010 (S)-VAL / L-VAL, with DNAn

Dimensions: 250 x 4 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Methanol

Pct. A: 60

Eluent B: Water

Pct. B: 40

Flow Rate: 0.6 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTE

1 1,1'-Bi-2-naphthol

2 1,1'-Bi-2-naphthol

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14143/CHX/J-03