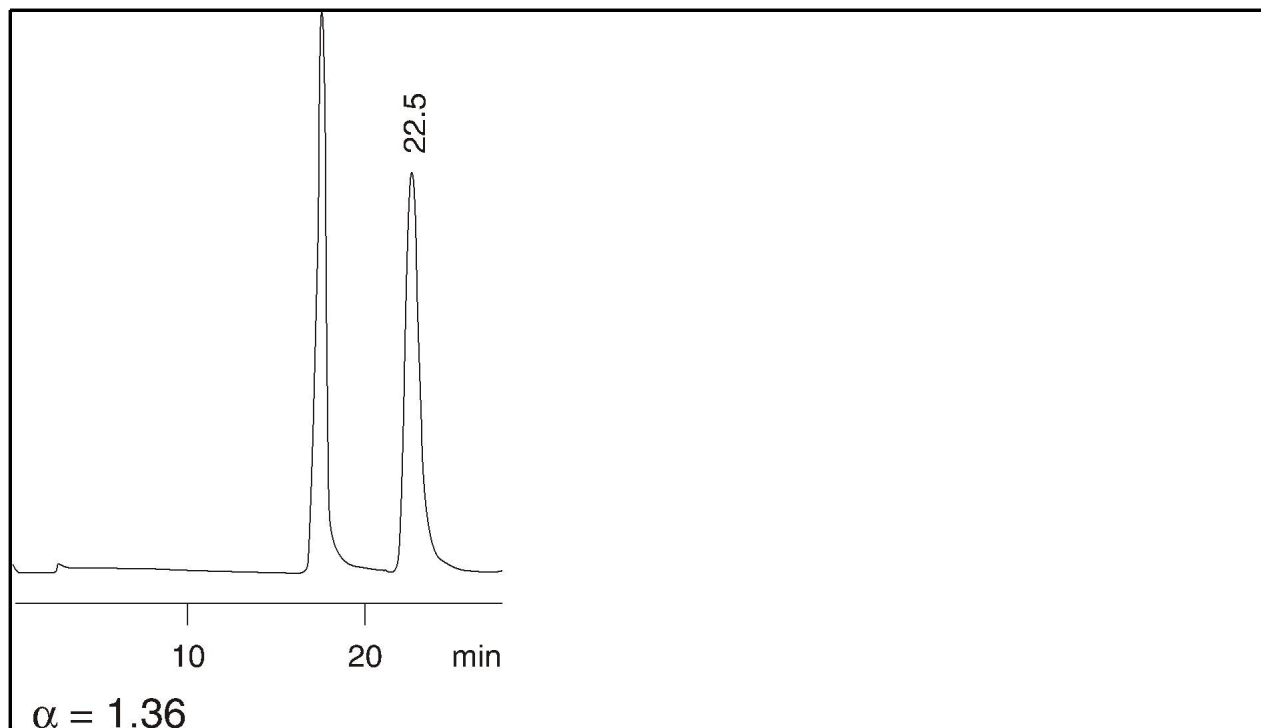


1,1'-Bi-2-naphthol

Column: Chirex 5 μ m Chiral I (Pirkle) 3020 (S)-t-LEU / L-LEU, with (R)-NEA

Dimensions: 250 x 4 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Hexane

Pct. A: 90

Eluent B: 1,2-Dichloroethane

Pct. B: 9

Eluent C: Ethanol

Pct. C: 1

Flow Rate: 1 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTE

1 1,1'-Bi-2-naphthol

2 1,1'-Bi-2-naphthol

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14144/CHX/J-03