


Benzalkonium chloride - Modified Method

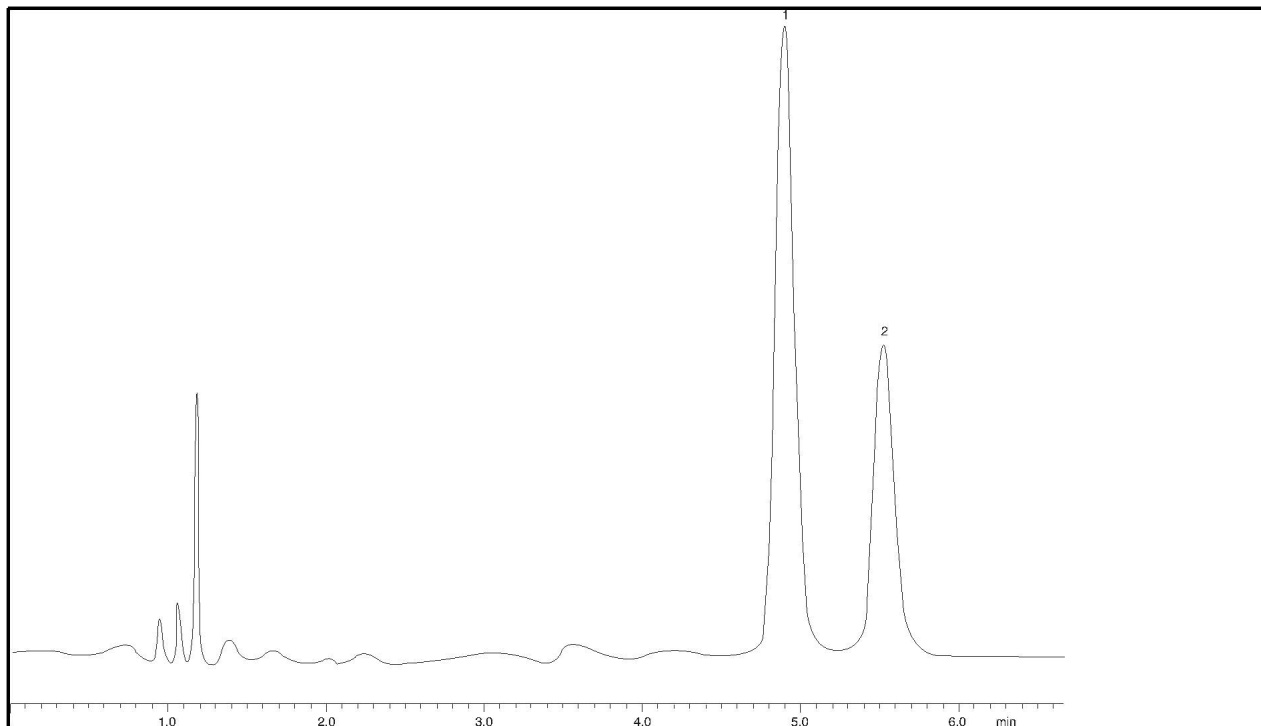

Column: Luna 5 μ m CN

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 100mM Sodium acetate pH 5.0

Pct. A: 40

Eluent B: Acetonitrile

Pct. B: 60

Flow Rate: 2 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTE

1 Benzalkonium chloride - C12

2 Benzalkonium chloride - C14

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 9300/USP/20