


Benzalkonium Chloride

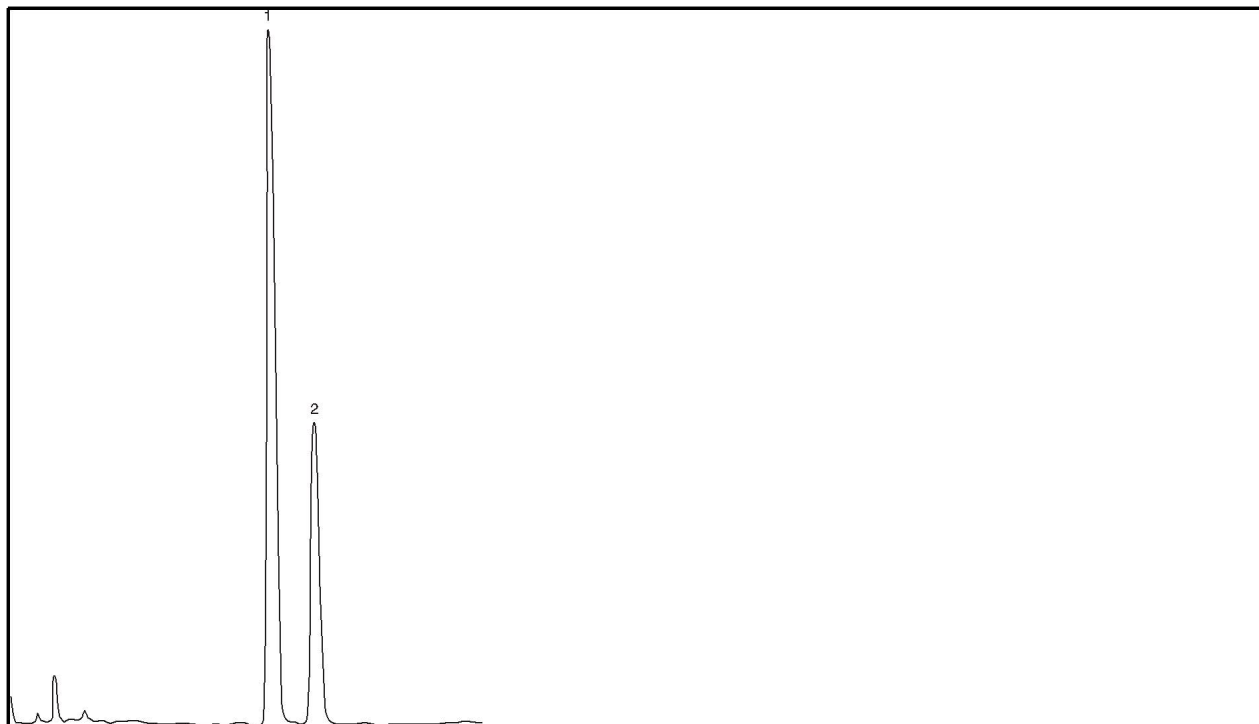

Column: Luna 5 μ m CN

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Acetonitrile

Pct. A: 55

Eluent B: Tetramethyl ammonium chloride buffer 50 mM, pH 3.0

Pct. B: 45

Flow Rate: 1 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 214 nm

SAMPLE ANALYTE

1 Benzalkonium chloride - C12

2 Benzalkonium chloride - C14

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 9905/ALG/202