


Estrogens on Luna Amino

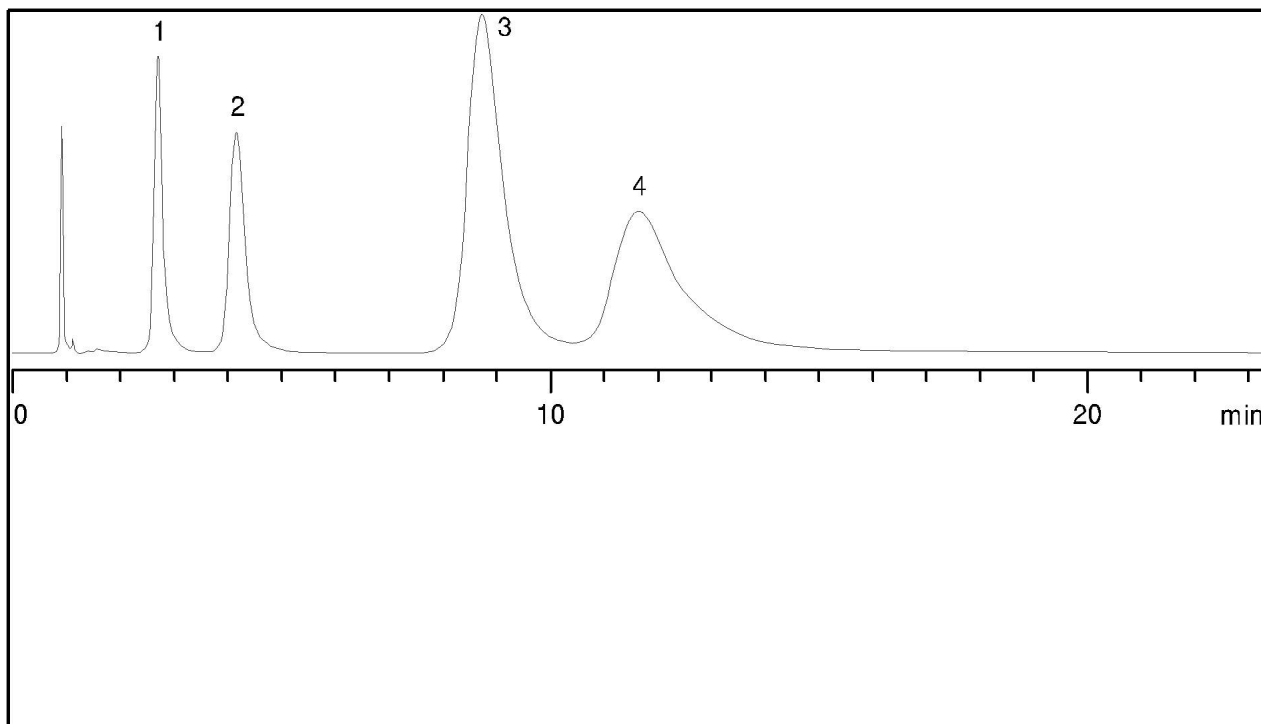

Column: Luna 5 μ m NH₂

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Hexane

Pct. A: 70

Eluent B: Isopropyl Alcohol

Pct. B: 30

Flow Rate: 2 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 280 nm

SAMPLE ANALYTE

- 1 Estrone
- 2 beta-Estradiol
- 3 Diethylstilbestrol
- 4 Estriol

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14298