


Ginger Root (*Zingiber officinale*)

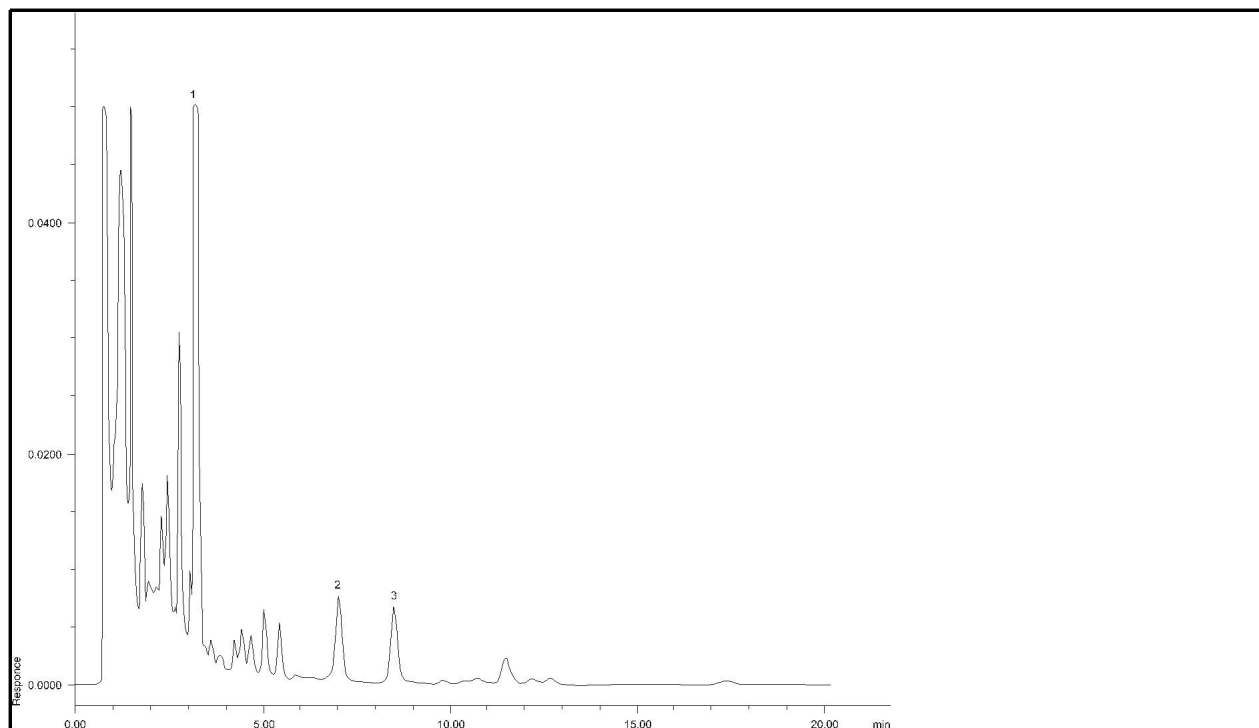

Column: Luna 5 μ m C18 / ODS (2)

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Acetonitrile

Pct. A: 50

Eluent B: Water

Pct. B: 45

Eluent C: Tetrahydrofuran

Pct. C: 5

Flow Rate: 1.5 mL/min

Temperature: 25 °C

Detection: UV-Vis Abs. - Variable Wave. @ 282 nm

SAMPLE ANALYTE

1 6-Gingerol

2 8-Gingerol

3 6-Shogaol

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1078/NUT/ginger2.doc