


## Proteins on Jupiter C4 in 10mM TFA

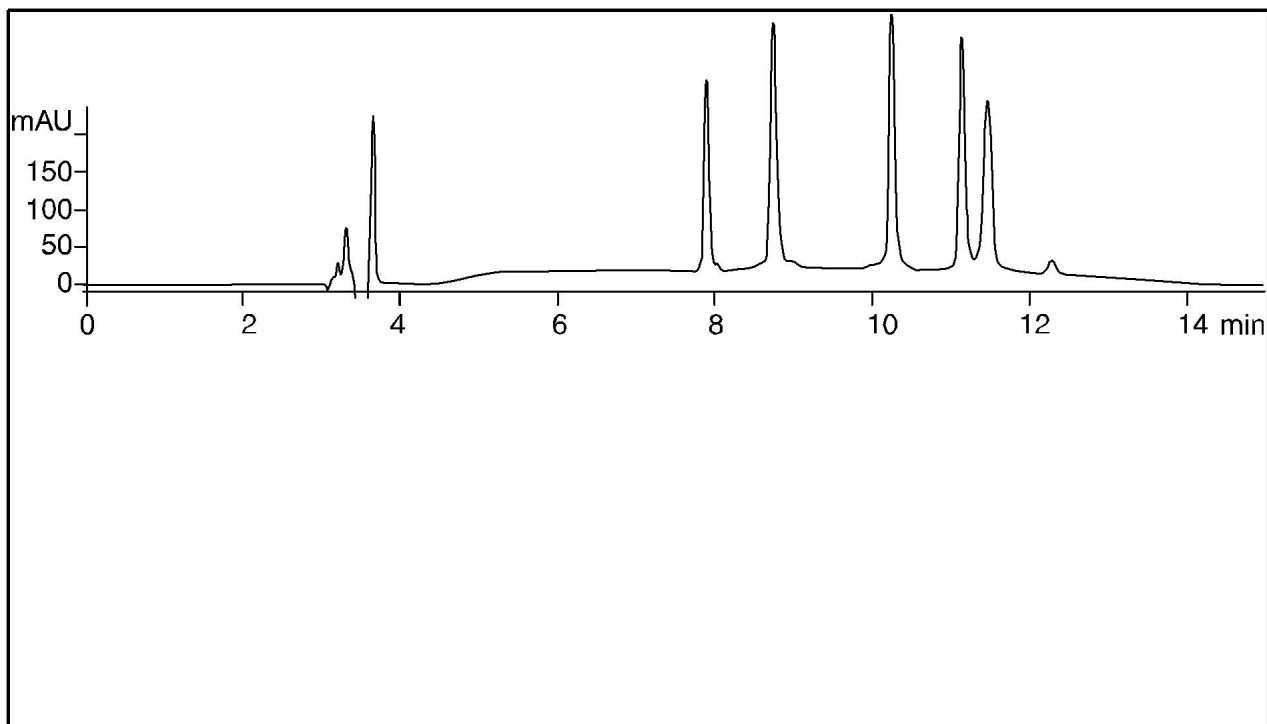

**Column:** Jupiter 5  $\mu$ m C4 300Å  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:**  
**Elution Type:** Gradient  
**Eluent A:** Water with 10mM TFA  
**Eluent B:** Acetonitrile with 10mM TFA

**Gradient Profile:**

| Step No. | Step Time | Pct A | Pct B |
|----------|-----------|-------|-------|
| 1        | 0         | 75    | 25    |
| 2        | 20        | 5     | 95    |

**Flow Rate:** 1 mL/min  
**Temperature:** 35 °C  
**Detection:** UV-Vis Abs. - Variable Wave. @ 214 nm

**SPE:****SAMPLE ANALYTE**

- 1 Insulin (12 kD)
- 2 Trypsinogen (24 kD)
- 3 Lactalbumin (14 kD)
- 4 Myoglobin (17 kD)
- 5 Carbonic anhydrase (29 kD)

**For more information call your Phenomenex Representative at (310) 212-0555**

Copyright 2000 Phenomenex Inc.

**Ref. No.:** 14184