


Mannitol and Sorbitol

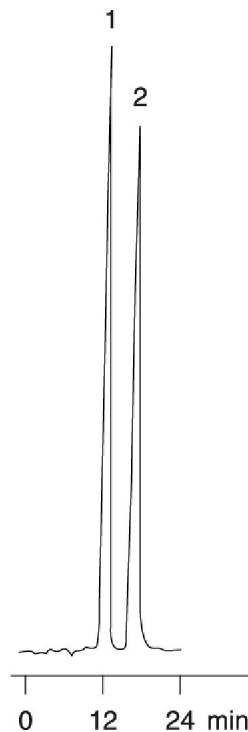

Column: Rezex 8 μ m Sulfonic Acid (SCX-SO₃H), strong cation exchange RCU-USP Sugar Alcohols

Dimensions: 250 x 4 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Water

Pct. A: 100

Flow Rate: 0.2 mL/min

Temperature: 85 °C

Detection: Refractive Index (RI) @ nm

SAMPLE ANALYTE

1 Mannitol

2 Sorbitol

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

©Copyright 2000-2004 Phenomenex Inc.

Ref. No.: 5516/CAT99/238