


Tryptic map of Dog Cytochrome C - Jupiter Proteo

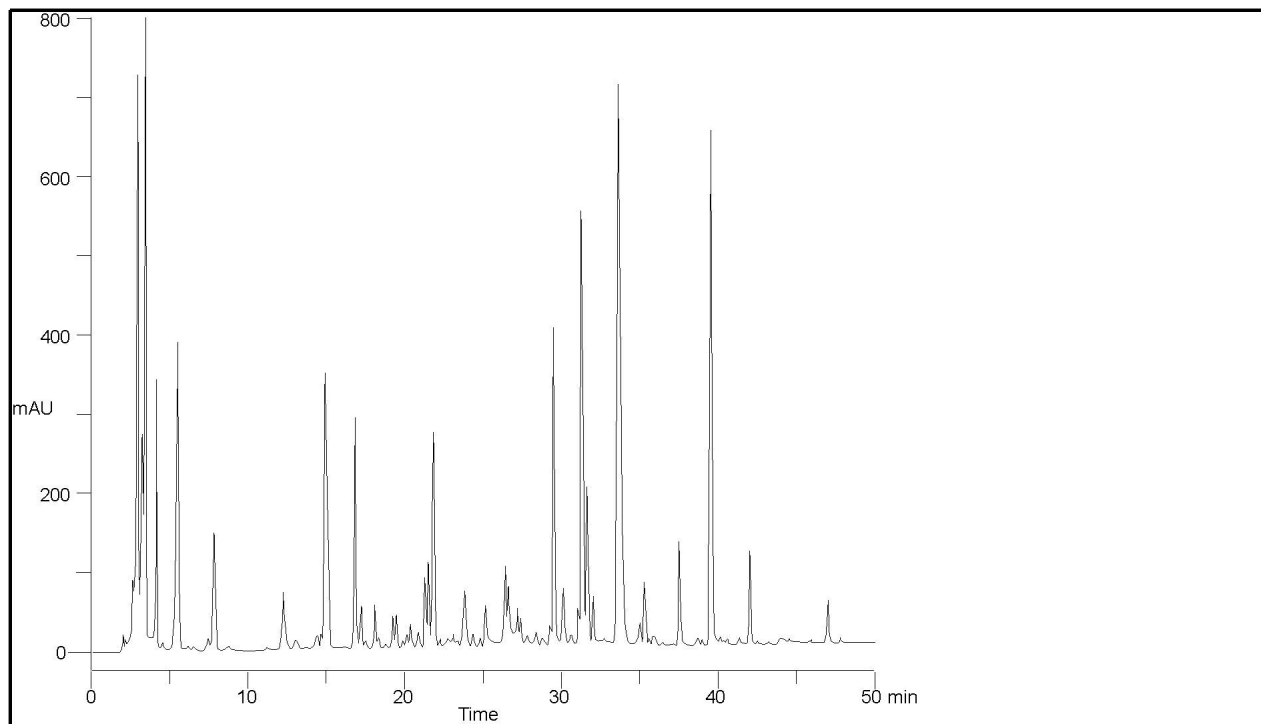

Column: Jupiter 4 μ m C12 Proteo
Dimensions: 250 x 4.6 mm ID
Order No:
Elution Type: Gradient
Eluent A: Water with 1.25mM TFA
Eluent B: Acetonitrile with 1mM TFA

Gradient Profile:		Step No.	Step Time	Pct A	Pct B
		1	0	95	5
		2	5	95	5
		3	60	60	40

Flow Rate: 1 mL/min
Temperature: 22 °C
Detection: UV-Vis Abs. - Variable Wave. @ 210 nm

SPE:

SAMPLE ANALYTE

- 1 Tryptic map of dog cytochrome c

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14370